



QTM

Qualisys Track Manager

GETTING STARTED

QUALISYS AB

Packhusgatan 6 · 411 13 Gothenburg · SWEDEN
Tel. +46 31 336 94 00 · Fax. +46 31 336 94 20
e-mail: sales@qualisys.com · www.qualisys.com



THE FIVE PICTURES ON THIS PAGE SHOW THE QUALISYS TRACK MANAGER SYSTEM IN ACTION. A VIDEO CONFERENCING, THREE-DIMENSIONAL IMAGE PROCESSING AND MOTION CAPTURE SYSTEMS IN ONE. THE QUALISYS TRACK MANAGER SYSTEM WILL FORM A UNIQUE LANDMARK IN MOTION CAPTURE TECHNOLOGY.

SPRING: "THE INFLUENCE OF THE SPRINGS ON THE BACK IN RIDING HORSES" (TELEVISION SPOT)

Qualisys AB, 2012-02-24

www.qualisys.se

No part of this publication may be reproduced, transmitted, transcribed, stored in a retrieval system, or translated into any language in any form by any means without the written permission of Qualisys AB

Copyright © 2011 Qualisys AB ® All Rights Reserved.

Qualisys Track Manager and QTM are registered ® trademarks of Qualisys AB

Table of Contents

QTM - Getting started	1
Introduction	1
System setup	3
Install QTM	3
Projects	3
Using projects	3
Projects with qualified users	4
Projects with students	4
Using projects on multiple computers	5
QTM startup procedure	5
Creating a new project	6
Switch project	7
Project view	8
Project data tree	8
Details	10
Backup of project settings	11
Project presets	11
Project folder	12
Setting up the Oqus system	13
Oqus connectors	14
Oqus display	16
QDS	17
QDS menu	17
Oqus network configuration wizard	18
Setting the aperture and focus	19
Connecting other hardware	20
QTM user interface	21
QTM user interface and windows	21
RT(Preview)/File mode	22
2D view window and features	22

3D view window and features.....	24
Trajectories in 3D views.....	26
6DOF bodies in 3D views.....	27
Trajectory info windows and menu.....	28
Data in Trajectory info windows.....	29
Select and move trajectories.....	30
Data info window and menu.....	30
Plot windows.....	31
Project options dialog.....	32
QTM capture settings.....	33
Capture settings overview.....	33
Connecting to the system.....	33
Loading linearization files.....	34
Calibration settings.....	34
Oqus marker settings.....	35
Calibration.....	37
Calibration overview.....	37
Camera positioning.....	37
How to perform a Wand calibration.....	38
Capture.....	41
Capture overview.....	41
Real time mode.....	41
How to use real time.....	41
Capture measurement files.....	43
How to capture a file.....	43
How to capture an Oqus video file.....	45
Batch capture.....	46
Processing.....	47
Processing overview.....	47
Batch processing.....	47
Identification of trajectories.....	49
Manual identification of trajectories.....	49
Automatic Identification of Markers (AIM).....	50
Generating an AIM model.....	50
Guidelines for data added to AIM models.....	53
Applying an AIM model.....	54
How AIM identifies trajectories.....	56
Data export to other applications.....	56
IR radiation notice and environmental regulations.....	59

IR radiation notice	59
EU customer information about about Waste Electrical and Electronic .. Equipment (WEEE)	59
Qualisys WEEE Collection Process	59
China ROHS (in Chinese)	61
有害物质声明	61

QTM - Getting started

Introduction

The Getting started manual includes the basics of how to use QTM with an Oqus camera system. For a detailed description of all of the features please refer to the **QTM User manual**. The complete manual can be found in PDF format on the entry for QTM in the Windows Start menu. It can also be found in HTML format in QTM on the Help menu.

On the Client login site of www.qualisys.com, you can find the current version of QTM, manuals, plug-ins and other information. Use your QTM user name and license key to log in. If you have any other questions please send it to support@qualisys.se.

The following subjects are described in the Getting started manual.

System setup

How to connect the Oqus system and set up the computer.

QTM user interface

Basic information about the QTM user interface, including a short description of the different windows and menus.

QTM capture settings

Information about the basic settings needed to be able to capture a measurement, including calibration setting and camera settings.

Calibration

How to perform a calibration.

Capture

How to make a capture/batch capture.

Processing

A short description of the different processing steps and a more detailed description of AIM (Automatic Identification of Markers).

System setup

Install QTM

Make sure you are logged in with an administrator account before you start installing QTM. To install the software, insert the QTM installation CD in the CD drive. If the installation program does not auto-start, locate and execute the setup.exe file on the installation CD. Follow the instructions given during the installation.

 Note: QTM will also install the DHCP server called QDS which is needed for the Oqus cameras

The first time QTM is started you must enter a user name and a license key. This is provided on the front cover of the QTM installation CD.



 Note: If the license is time limited you must check the **Time limited** checkbox and enter the correct expiration date.

If there is an internet connection, QTM will automatically check for updates when it is started. You can also use the **Check for updates** option on the **Help** menu. If the QTM computer does not have an internet connection, you can look for the latest software from another computer on the web page www.qualisys.com. Log in with your user name and license key. The user name and license key can also be found in QTM. Go to **About Qualisys Track Manager** on the **Help** menu.

Projects

QTM needs a project to capture measurement data. The project is a folder that contains all the files and information needed for QTM to process the data. A project can therefore be easily transferred to for example another computer with all the settings and files needed for the processing. To create and use projects follow the instructions in the chapters below.

Using projects

The projects can be used in a lot of different ways. The following is advice on how you can use the projects, but it is recommended to read about the different and adapt to your needs.

Think about the following when deciding how to use the projects.

1. When to create a new project.
 - a. It is usually recommended to use one project per marker setup. Then you can save different subject in sub-folders of the project data folder.

- b. You may of course create any amount of projects based on for example the people using the lab.
2. Settings to use when creating a new project, see chapter “Creating a new project” on page 6 and “Project presets” on page 11.
3. It is recommended to save all of the measurement in the project data folder, because then you can browse them in the Project view window, see chapter “Project view” on page 8.
4. Do you want to backup settings, see chapter “Backup of project settings” on page 11.
5. Do you want to view the list of recent projects at startup or always open a certain project, see chapter “QTM startup procedure” on the facing page.

Read the following sections for advice on specific setups: “Projects with qualified users” below, “Projects with students” below and “Using projects on multiple computers” on the facing page.

Projects with qualified users

If there is more or less one user or multiple qualified users of QTM and the system then the setup of the projects can be done as the following example.

1. Create a new project when there is a new study, with a new marker setup.
2. Save all of the QTM files in the project data folder as you can then browse the data in the Project view in QTM, see chapter “Project view” on page 8.
3. Make a backup of the project settings if you like to be sure that you can always go back to the settings you know are correct, see chapter “Backup of project settings” on page 11. If someone opens the project any changes they do to the settings are saved automatically.
 - a. If you want to completely protect your projects from other users it must be done via Windows, for example by having multiple logins to the computer.
4. Pin projects that are often used in the **Open project** dialog, see chapter “Switch project” on page 7.
5. Create a **Project preset** if you create a lot of new projects and want to use the same settings when it is created, see chapter “Project presets” on page 11.

Projects with students

If there are a lot of students learning QTM then the setup of the projects can be done as the following example.

1. Create a Project preset so that the students can start with the same settings, see chapter “Project presets” on page 11.
2. Tell the students to create a new project with the Project preset and save them at a specific place. The default folder is Documents in Windows, but it can be changed, see chapter “Folder options” in the QTM user manual.
 - a. It is a good idea to use a different Windows login for the students so that they cannot access other projects.

- b. If the students can access the other projects make sure that you backup the settings of the projects, see chapter “Backup of project settings” on page 11.
- 3. Make sure that the QTM files are saved in the project data folder. It is the default path if nothing has been changed in the project settings.

Using projects on multiple computers


If you use a project with different camera systems you can either move the whole project (or computer) or just move the Data files and AIM models.

- It is recommended to setup the project on both computers and just move the Data and AIM files. Because then there is no risk for losing settings, on the other hand you have to remember to change the settings on both places if you need to.
 - If you move the whole project, then you need to decide on how to handle the camera settings. QTM will automatically detect on **New** if the cameras have changed and you will get the option to import all the settings from another project. It is important to notice that it will overwrite all the settings in the current project, so you will need to add any settings that differs again, for example the AIM model.

 Note: If you locate the system first, QTM will not give a warning for a new camera system, because it is seen as you then wanted to check for a new system.
 - If you just move the Data files there is no problem in opening them on the other computer, because all of the information needed to reprocess the file is included in the file.

QTM startup procedure

The first time you start QTM on a new computer you will get a question about creating a project.

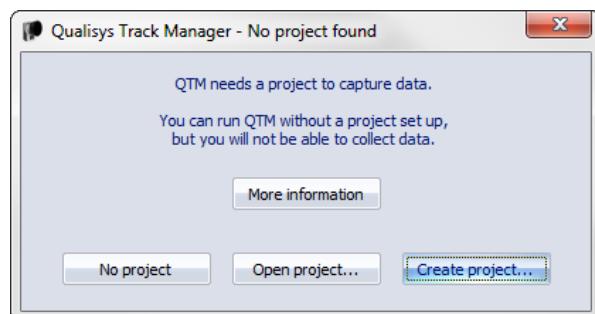

There are three ways to start QTM when you have not used projects before.

Create project

This is the default option because you must have a project to capture data in QTM, see chapter “Creating a new project” on next page.


Open project...

Use this option to open a project folder that has been copied from another computer.

No project

If you only want to open QTM files you can start QTM without a project, but you will not be able to capture any data.

After using projects the first time the default startup procedure for QTM is to display the **Open project** dialog. Each time you start QTM you have to select the project to load, to make it less likely that someone changes the projects settings by mistake. For more information about the dialog see chapter “Switch project” on the facing page.

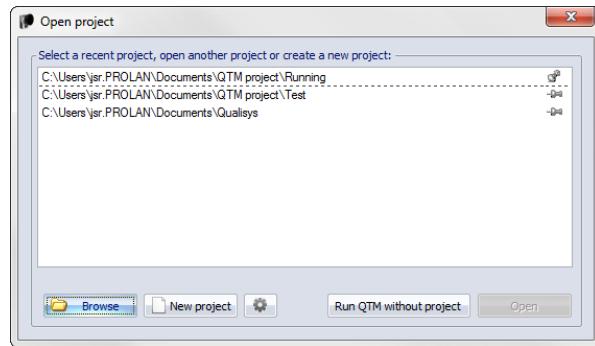


It is also possible to choose to open the most recent project or a selected project

when starting QTM. To set these click the button, see chapter “Startup” in the QTM user manual.

Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

Creating a new project

A new project is created with the **New project** option on the **File** menu or in the **Switch project** dialog. The following dialog is opened with options to create a new project.

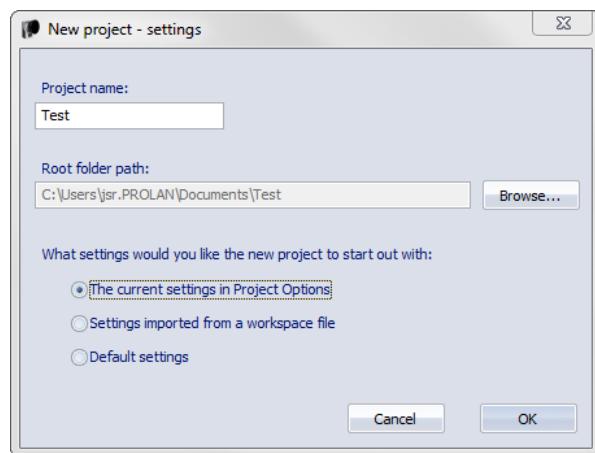

Enter the following information to create the project

Project name

The project name will also be the name of the project folder . The folder is created automatically if it does not exist at path specified in **Root folder path**.

Note: To change the project name after it has been created just change the folder name.

Root folder path

This is the path of the project folder, the default is the Documents folder in Windows. Click **Browse** to change the path.

 Note: The **Default path for new projects** can be set on the **Folder options** page in the **Project options** dialog.

 Note: If you want several computer users to access the same project, it must be saved in a folder which all users have access to. For example on Windows 7 you can use the folders of All users.

Base the new project on:

Select the source of the settings that are used when creating the project with the drop-down box.

Custom preset: ...

Any project presets available on the computer are displayed at the top of the list, for more information about the presets see chapter “Project presets” on page 11.

 Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

The current settings

The settings are copied from current settings in **Project options**. This is the default option.

Settings imported from another project

The settings are imported from another project. After clicking OK, you have to browse to the project settings file that you want to use. Any settings that is not set in the workspace file will be set to default values.

 Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

Settings imported from a workspace file

The settings are imported from a workspace file saved in QTM 2.5 or earlier. Any settings that is not set in the workspace file will be set to default values.

 Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

Default settings

All settings are set to the default values and all of the current settings are deleted.

IMPORTANT: If you select any of the options other than the **Current settings**, you will lose the current settings in **Project options**. If you haven't got any projects in the recent projects list, then the settings are saved in a backup in c:\ProgramData\Workspace backups (Windows 7) or c:\Documents and Settings\All Users\Application Data\Qualisys\Workspace backups (Windows XP).


Switch project

Use **Switch project** on the **File** menu to manage the projects, for example to switch to a recently used project or open another project.

 Note: Opening a project only loads the settings in **Project options** and other settings such as those in the **Start capture** dialog. QTM files outside the project can still be opened and processed if you like. QTM will then note that the file is not in the current project in the title bar.

 Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

The calibration is also checked if the cameras are changed at **Locate** or **New**. If there is a matching calibration for the new cameras it will be loaded in the project.

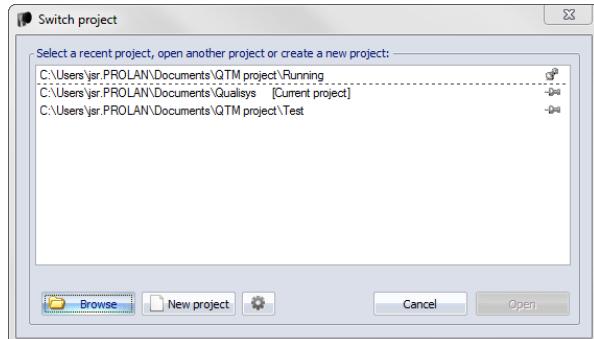

The list in the **Switch project** dialog displays the 100 most recently used projects. Double-click on a project to open it. You can pin  projects so that they are always on the top of the list, you have to click **Open** for the pinning to change.

The current project is shown with **[Current project]** next to the name. The current project is also displayed in the title bar of QTM and **Project options**.

 Note: When QTM starts up it checks if the projects exist on the computer and purges the list of any deleted projects.

The following options are available in the dialog:

Browse

Use this option to open a project that is not in the list above, e.g. one copied from another computer.

New project

Create a new project, see chapter “Creating a new project” on page 6.



Open settings for the startup of QTM, see chapter “Startup” in the QTM user manual.

Open

Open the project selected in the list above.


Project view

The **Project view** window displays the data files included in the current project. The view consists of two parts: **Project data tree** and **Details**. The name of the current project is displayed at the top of the Project view.

To open the **Project view** window go to **View** menu and select **Project view**, the shortcut to toggle the window is **ctrl+r**.

Project data tree

The **Project data tree** window is used to display and open the files in the Data folder of the current project. This includes folders and QTM files, but also all other files.

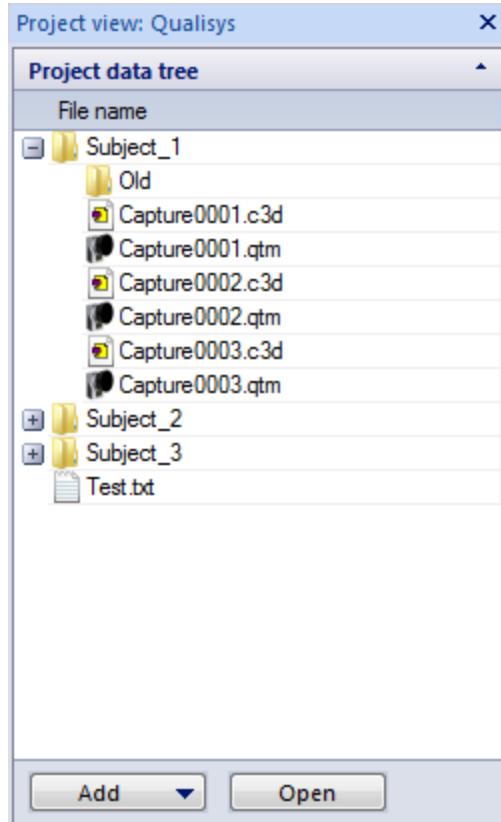


There are two commands at the bottom of the data tree: **Add (New folder or New measurement)** and **Open**.

When adding a folder it is placed in the currently selected folder, if none is selected it is placed in the root of the Data folder. Adding a new measurement is the same as starting a new capture with the **Capture** command on the **Capture** menu and the folder to save the file in is set to the currently selected folder in the project view.

The **Open** command is only available when a file is selected and it will then open that file. If it is a QTM file it is opened in the current QTM window or a new window, depending on the current settings on the **GUI** page in **Project options**. Other files are opened with its corresponding Windows program.

The QTM file that is open in the current QTM window is displayed in bold. If you open a file that is already open in another QTM window you will switch to that window, unless you use **Open in New window** and then you will open a new instance of the file.

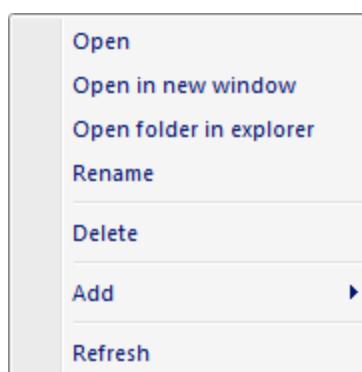

Right-click on a file or a folder to get the following options.

Open (only available for files)

Opens the selected file in QTM if it is a QTM file. Otherwise it is opened by its corresponding Windows program.

Open in New window (only available for QTM files)

Opens the QTM file in a new QTM window. It overrides the **Close current file** option on the **GUI** page.

Open folder in explorer

Open the currently selected folder in Windows explorer so that you can for example move the files.

Rename

Rename the selected file or folder.

Delete

Delete the selected file or folder.

Add

New folder

Add a new folder in currently selected folder. If no folder is selected it is placed in the root of the Data folder.

New measurement

Add a new measurement, which is the same as **Capture** on the **Capture** menu and the folder to save the file in is set to the currently selected folder in the project view.

Refresh

Refresh the content in the **Project data tree** so that it matches the content of the project data folder. This is usually done automatically by QTM.

Details

Details for item 'Capture0001.qtm':	
Property	Value
File name	Capture0001.qtm
File type	QTM motion capture...
File size	1815 KB
File created	2011-11-21 21:06
File modified	2012-02-01 13:15
Full path	C:\Users\jsr.PROLA...

The **Details** window displays information about the currently selected file. The displayed data is:

File name

File type

File size

File created


This is not the time when the data was captured, but when the file was created on the computer which is not the same if the file has been copied.

File modified

Full path

Backup of project settings

The project settings can be saved in a backup file. The backup can then be used if you want to use two different sets of settings in the same project or if you want to be sure that the settings are not changed by someone else, because the project settings are saved automatically when they are changed in QTM. Click on **Settings backup and restore** in the **File** menu to create and manage the backups.

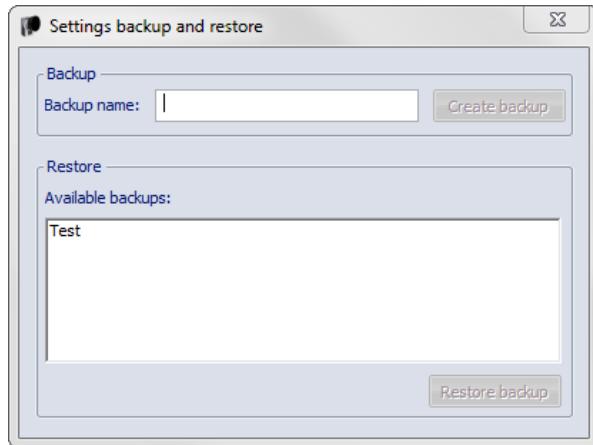

Backup

Enter the **Backup name** and click **Create backup** to save a backup of the current settings in **Project options**.

Restore

The list displays all of the backups saved in the current project. Select one and click **Restore backup** to copy the settings to **Project options**.


 Note: The current settings in Project options will be replaced.

 Note: QTM will use the latest calibration made on the computer, that was made with the same cameras (placed in the same order). It means that the calibration may be loaded from another project when switching projects, it is then copied to the Calibrations folder in the project folder.

Project presets

The project presets can be used when creating projects to make sure that you start with the same settings. The project preset contains all of the settings QTM (Project options and other settings such as those in the Start capture dialog). The preset also contains any AIM model that was in the project AIM folder when creating the preset.

 Note: The presets are saved globally so that it can always be accessed when creating projects in QTM.

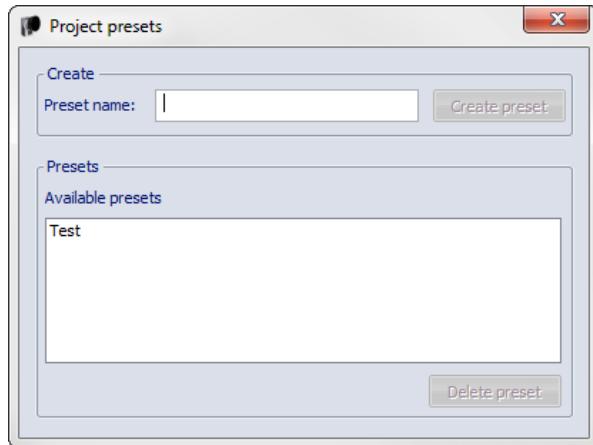

To create a preset follow these steps.

1. Open a project with the settings and AIM models that you want to use.
2. Go to the **File** menu and select **Settings management/Project presets** to open the **Project presets** dialog.
3. Enter the name of the preset and click on **Create preset**.

To use a preset follow these steps.

1. Go to the **File** menu and select **New project**.
2. Enter a **Project name** and select the settings to **Base the new project on** from the drop-down list. The presets are listed as **Custom preset:** followed by the name.

A preset can be deleted from the **Project presets** dialog. It is important to notice that the presets are deleted from the computer so that it can no longer be accessed by anyone.

Project folder

The project folder is used for organizing a QTM project. The folder is displayed in Windows explorer with a special icon to indicate that it is a QTM project folder. It contains the following files and folders:

Data

This is the default location for the captured QTM files. You can create subfolders in this folder if you want to sort the files, for example for different subjects.

AIM models

This folder contains all of the AIM models created in the current project.

Calibrations

This folder contains all of the calibrations made in the current project.

Labels

This is the default folder for Label lists.

Settings

This folder contains the backups of project settings.

Settings.qps

This file contains the current settings of the project.

QTM creates a default project folder called Qualisys in My Documents when installing QTM on a new computer. However you can create a new project folder anywhere. For example if you want different users on the computer to access the

same project it can be saved under the Public user in Windows 7. Copy the whole folder if you want to share the project with someone else.

There are other settings and files which are needed for all users that use QTM on a computer. To see where these are saved check the **Folder options** page in the **Project options**.

 Note: The QTM program and other components installed by the QTM installer are placed in Qualisys Track Manager folder under \Program Files\Qualisys (Program Files (x86) for Windows 7 64 bit).




Setting up the Oqus system

Before you connect the Oqus camera system, make sure that the QDS (Qualisys DHCP server) is running and that the network interface settings are correct. For more information about QDS see section “QDS” on page 17.

The Oqus system is easy to setup. The connectors are unique and cannot be connected to the wrong ports. Further, the connector color matches that of the port. The **DATA** connector can be connected to any of the two **DATA** ports, and the **POWER** connector can be connected to any of the two **POWER** ports, so it does not matter on which side you put the connector. For more information on the connectors, see section “Oqus connectors” on next page.

 Note: When the cables have been connected correctly the LEDs on the back of the Oqus will be lit. The **EXT** LED will be lit green and the **ACT** LEDs will be blinking.

However, the connection of the power adaptors does require some attention. One power adapter can power up to 4 Oqus cameras. Therefore, the connection of a camera system comprising of more than four cameras must look something like the image below.

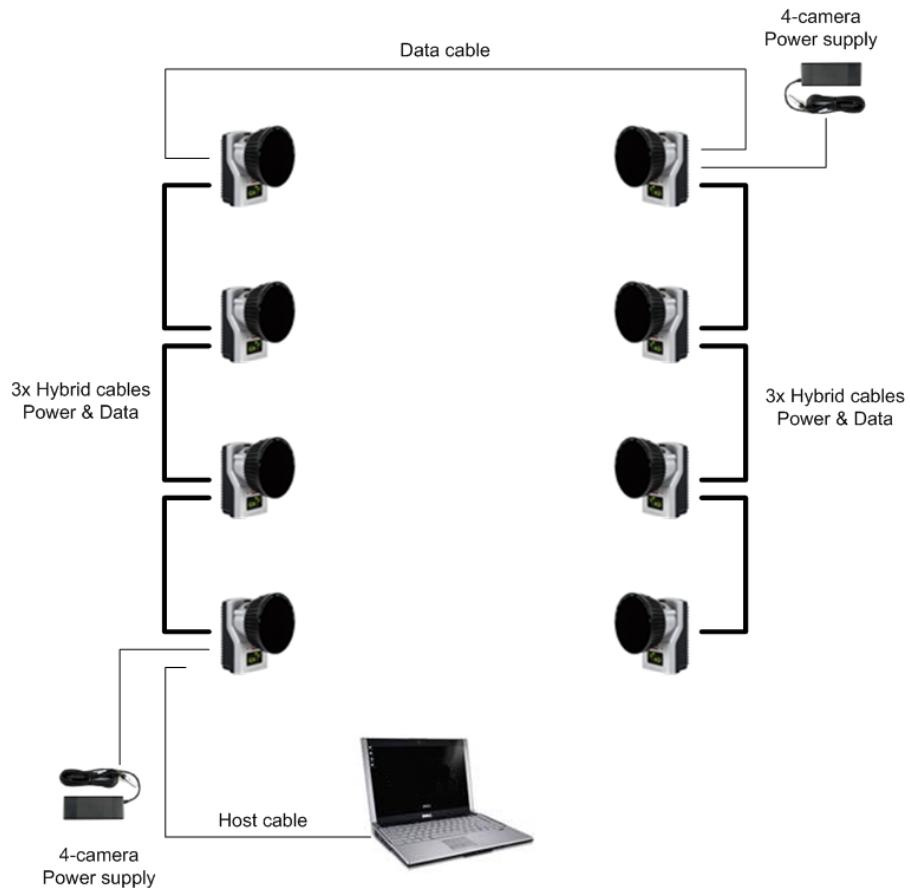


This means that you must use the following cables for an 8 camera system:

- 1 host cable between computer and camera.
- 6 bundled cables with power and data.
- 1 data cable connected between camera 4 and 5 in the setup, because one AC/DC adapter can only power 4 cameras.
- 2 AC/DC adapters, connected to for example camera 1 and 5.

Oqus connectors

The back of the camera holds six connectors for power, data and control connections. The view differs slightly depending on the type of camera. The image below shows the standard version of the camera. The water protected version uses different connectors and lacks the LEDs found on the standard version.

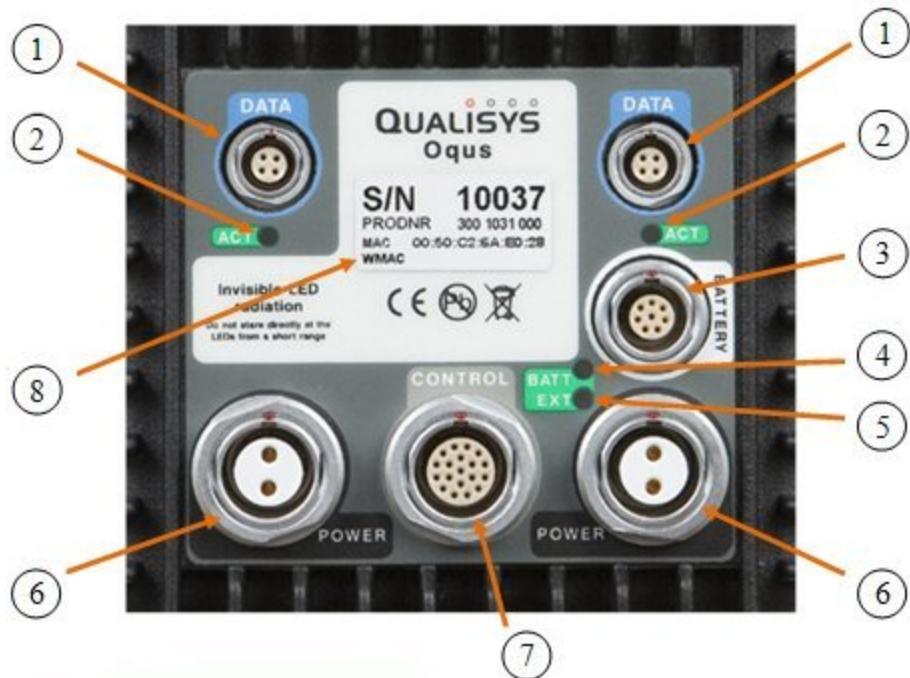


1. **Left Data port, Right Data port** (light blue)
Ethernet connector. 100BaseTX/802.3i, 100 Mbps, Fast Ethernet. The two ports are identical.
2. **Left Ethernet activity indicator, Right Ethernet activity indicator**
Shows the status of the Ethernet connection. Fixed green light means that a carrier signal has been detected and that the connection is up. Flashing green light indicates that data is received and/or transmitted.
3. **Battery port** (white)
Used to supply the camera with power from an Oqus compatible battery.
4. **Battery status indicator**
Lit green when the camera is supplied through the **BATTERY** port.
Lit red when a voltage outside the specified range (10-16V) is connected to the port.
5. **Power supply status**
Lit green when the camera is powered through one of the **POWER** ports.
A red light indicates internal power supply error.
6. **Left power supply port, Right power supply port** (black)
Daisy-chain power port. Supplies the camera with 48VDC and can be daisy chained to supply cameras further down the chain with power. The ports are identical there is no specific in-port or out-port.
7. **Control port** (light grey)
The control port is used to synchronize the camera with external sources, and contains pins for e.g. external trigger in, external sync in and external sync out. Splitter cables are needed to connect one or more BNC cables to this port.
8. **Camera identification**
This label provides information on:
 - The serial number of the camera
 - The product number

- The Ethernet Mac address
- The WLAN Mac address

Oqus display

The Oqus camera has a large graphical OLED display and three LEDs on the front to inform the user of the current status of the camera. The display shows, among other things, the camera number and the number of markers currently seen by the camera.

 Note: The display will be turned off when the camera enters stand-by mode, i.e. if the camera has not been in use for 2 hours. Start a preview in QTM to light up the display again.

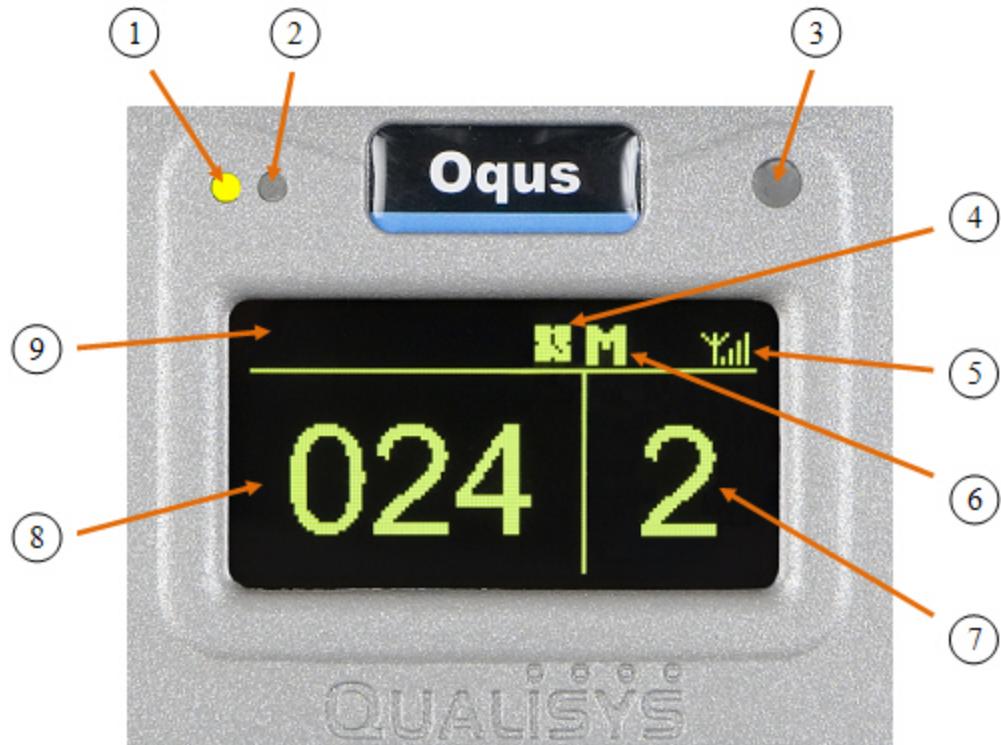

1. Measurement status indicator

- Green light - The camera is ready to start a measurement
- Yellow light - The camera is measuring
- Flashing green light - Waiting for trigger to start measurement
- Flashing yellow light - Waiting for trigger to switch from pre-trigger to post-trigger measurement

2. Error indicator

A red light indicates that an error has occurred. The LED is blinking when a software error occurs and is lit constantly if a hardware error occurs.

3. IR receiver

The IR receiver is used for synchronization with certain active markers. It detects modulated light with a frequency of 33 kHz and is sensitive to light with wavelengths between 800 and 1100nm.

4. Synchronization status

During the synchronization phase this symbol is flashing. When the camera is synchronized with the master camera in the system it becomes stable.

5. WLAN indicator

This symbol is displayed when the WLAN of the camera is activated.

6. Master/Slave indicator

An M indicates that the camera is master for the system and by that controls for example internal synchronization. An S indicates that the camera is a slave. The indicator can also be a rotating + sign, which means that the camera is looking for the Master camera.

7. Camera number

The area to the right usually shows the camera number that the camera has in QTM. The camera number can be changed with the **Reorder** tool in the 2D view window. This number is stored in the camera so it is shown at the next camera startup.

 Note: If the camera has never been connected to QTM the last three digits of the serial number (upper part) and the last octet of the IP-number assigned to the camera (lower part) will be shown instead. This can also be activated from the QDS menu, see "QDS" below.

8. Marker area

During a marker measurement this area shows the number of markers currently seen by the camera. When the camera is idle or is collecting video, this area shows '----'.

9. Text area

This area is used for scrolling text messages, for example during startup.

QDS

QTM comes with a DHCP-server called QDS (Qualisys DHCP Server), which distributes IP addresses to the Oqus cameras. An IP address for each camera is required to be able to communicate with them over the Ethernet network. QDS will be installed automatically with QTM and it must be running at all times, to provide the cameras with IP addresses at startup. The DHCP server will only give IP addresses to Oqus cameras so it will not disturb your computer network.

QDS includes a wizard to configure network interfaces to use them with Oqus cameras and to activate wireless operation on an Oqus camera.

QDS menu

To open the QDS menu, right-click on the QDS icon  in the status toolbar.



The following are the most important options on the QDS menu:

Oqus configuration wizard

The Oqus configuration wizard can be used to setup your network interface for Oqus and to configure an Oqus camera for wireless operation.

Advanced

Advanced configuration of network interfaces on the computer.


Network configurations

Using this sub-menu you can save or load network configurations. Click on **Save** to save the current configuration and then to load a configuration click on **Load**.


Oqus network configuration wizard

The Oqus network configuration wizard will guide you through the different steps to setup the network for Oqus cameras. Follow the steps below.

1. Click on **Oqus configuration wizard** in the QDS menu to start the wizard.

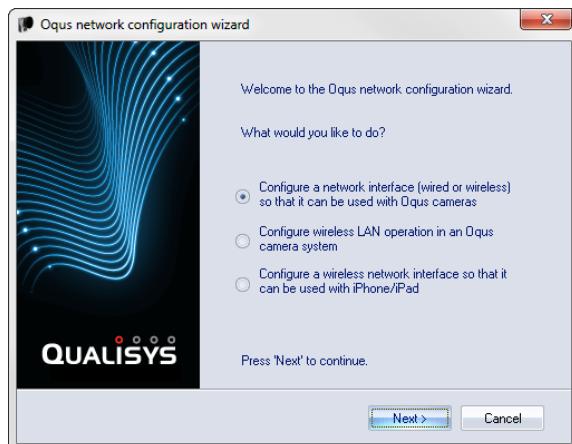


2. Select **Configure a network interface...** and click **Next**.

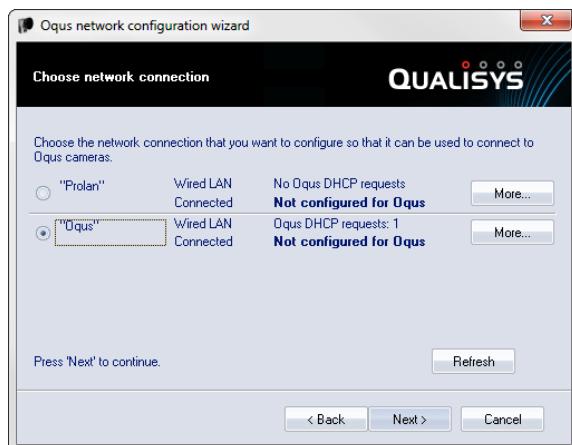


3. The list will show all enabled network interfaces on the computer. Select the interface that you want to use with Oqus cameras and click **Next**.

Note: The wizard will not configure all interfaces, for example a network interface that is already connected to an internal network and has received an IP address will not be configured because it is considered to have a running DHCP server. However any disconnected interfaces can be configured by the wizard.

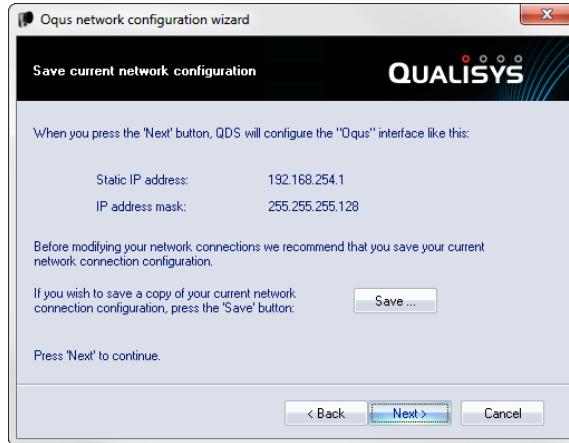

4. The wizard shows how it will change the selected interface. You can save the current network setup with the **Save** button for backup. Click **Next** to continue.

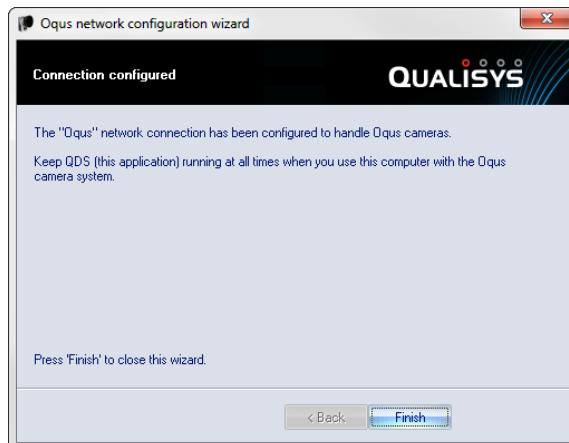

5. Click **Finish** to close the wizard.

Setting the aperture and focus

For Oqus it is very important to set the aperture and focus correctly for your measurement volume. If these are not set correctly the cameras will not reach their full potential in detecting markers. Follow this procedure to set the aperture and focus.



1. Turn the strobe part of the camera counterclockwise to expose the lens for adjustment.

2. Open the aperture as much as possible, i.e. set the lowest aperture.
 Note: For the Oqus 1 series and 3+-series it is sometimes possible to have the aperture on a higher number than the lowest. Check this by using the **Marker intensity** mode in QTM to see whether the markers are bright enough. The advantage with a more closed aperture is a longer focal depth.
3. Place a marker in the measurement volume. Use the same size as will be used in the actual measurement.
4. Switch the cameras to **Video** mode in QTM and change the focus until the markers are as small and sharp as possible. Make sure that the **Flash time** is long enough on the **Camera settings** page.
 - a. When you have set the focus correctly, switch to **Marker intensity** mode in QTM to make sure that the markers are visible in **Marker** mode. If the color of the markers is not red, set a longer exposure time in QTM on the **Camera settings** page in QTM.
5. Turn the strobe back (clockwise) to close it. The strobe should always be closed during measurement, to achieve the best strobe light distribution.


Connecting other hardware

Other hardware that can be connected to QTM includes analog boards, force plates and EMG. For description of how to do this and for the settings that need to be set, please refer to the [QTM User manual](#).

QTM user interface

QTM user interface and windows

The QTM user interface and windows are flexible and user-friendly. There are many features available in QTM and below you will find the main features described.

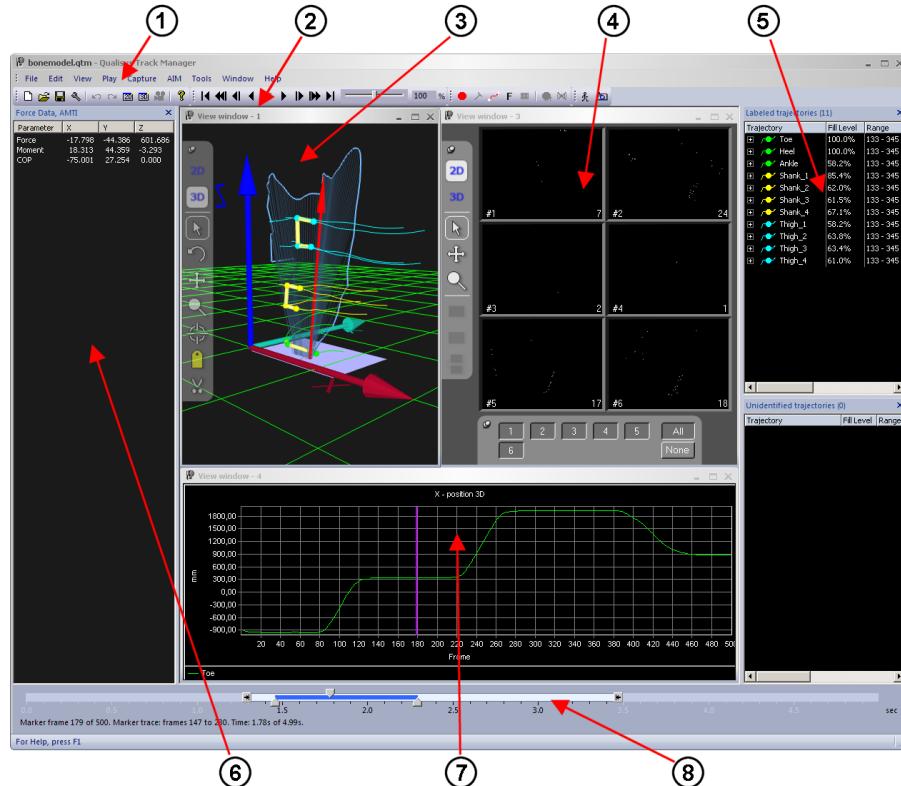

1. Menus

From the menus you can reach all of the features in QTM.

2. Toolbars

The toolbars contain buttons for the more frequently used features. The most important buttons are:

New file 

Project options 

Start capture 

Reprocess file 

Apply AIM model 

Generate AIM model 

3. 3D view window

The 3D view window displays the captured data in 3D. The data can also be edited by using the different tools in the **3D view** toolbar

4. 2D view window

The 2D view window displays the captured data in 2D for each camera. In RT/preview mode there is also a **Camera settings** sidebar where you

can modify the camera settings. There are also some tools in the **2D view** toolbar for modifying the camera in RT/preview mode. At the bottom of the window there are buttons to select which cameras that are displayed.

5. **Trajectory info windows**

The **Trajectory info** windows contain the tracked 3D data.

6. **Data info window**

The **Data info** window can display the following data: 2D, 6DOF, Analog and Force.

7. **Plot window**

Any of the data in the **Trajectory info** windows and the **Data info** window can be plotted.

8. **Timeline toolbar**

The **Timeline** toolbar controls the position in a file.

RT(Preview)/File mode

The QTM software has two modes of operation: RT (preview) mode and file mode.


The RT (preview) mode is entered when a new file is opened in QTM. Click on **New** on the **File** menu to open a new file. In RT (preview) mode the motion capture is shown in real-time before and during a measurement. In RT, the data can also be accessed via a real time TCP/IP protocol.

The file mode is used when the motion capture is finished. Manual editing of a QTM file can only be done in file mode.



2D view window and features

The 2D viewwindow display the raw data of all cameras connected to QTM. This includes 2D views of the motion capture data from each camera in the measurement system. It also includes Oqus video that are displayed in the 2D view of the Oqus camera and the DV/webcam video which is displayed after the motion capture cameras.

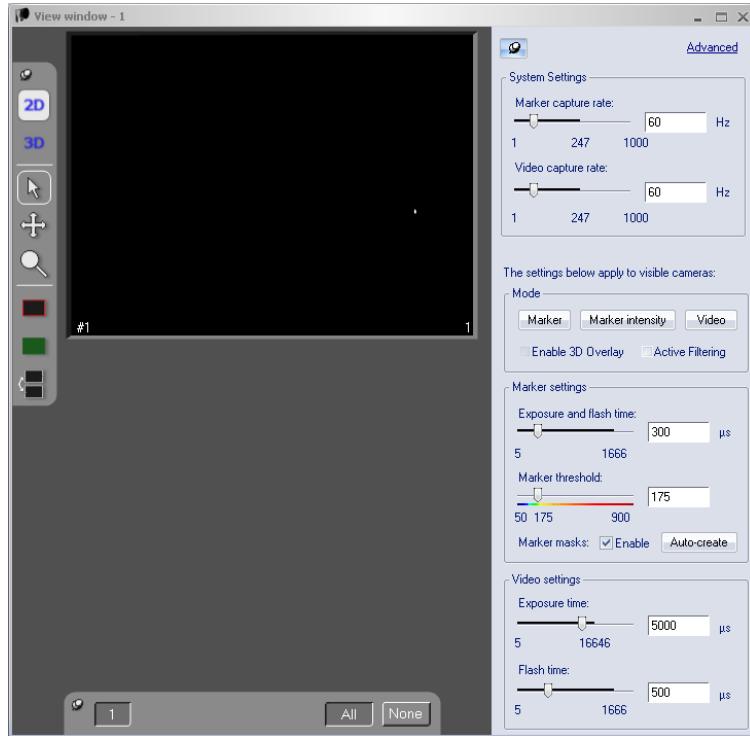

The 2D view window also includes the following controls:

2D view toolbar

Contains tools for the 2D view. The first three are for select, zoom and translate. Further, there are three tools for FOV, camera masking and camera order. Use these to change the image size, place marker masks to cover unwanted markers and for rearranging the order of the cameras.

Camera selection bar

Select the cameras that you want to view in the 2D view window. When changing the settings on the Camera settings sidebar it only affects the currently visible cameras.

Camera settings sidebar (Oqus)

Control the basic settings for Oqus cameras, for example you can change **Exposure time** and **Marker threshold**. Use these settings during the setup of the system, see section “Oqus marker settings” on page 35.

With the **Mode** options you select how to view the image from the camera during preview.

Markers

In this mode the cameras calculate the position of the markers.

Marker intensity

This mode shows a color coded image of the marker calculation and is useful when changing the camera settings.

Video

This mode shows the video image from the camera. If you want to capture video instead of markers from a certain camera, set it to video mode before starting the capture.

2D view window menu

Right-click a camera view to open the options for a single camera. Among other things you can rotate the 2D view to match the physical setup.


The following objects are displayed in a 2D view for a camera:

- The number in the lower left corner of the 2D view of a camera is the camera id. A motion capture camera is displayed as for example #1, and a DV/webcam camera is displayed as for example #1V.
- The number in the lower right corner is the current number of markers seen by the camera. A 'V' is displayed instead, if a motion capture camera is in video mode.
- The current **Image size** of the camera is shown as a red square and the part of the image that is outside the image size is greyed out. For Oqus cameras the **Image size** can be changed with the **Image size** tool on the **2D view** toolbar
- For Oqus the marker masks are shown as green squares.
- The Oqus video is also shown in the 2D view window, both in preview and in a file. This means that all of the actions, like zoom and 3D overlay, can be performed both in preview and in a file. Set an Oqus camera in Video mode and start a new capture to capture a video file from that camera.
- The DV/webcam video are displayed after the motion capture cameras in the 2D view window. You can use zoom on the video in both preview and file mode. The video cameras will appear in the same order as they are on the **Video devices** page in the **Project options** dialog.
- The 2D markers are color coded when marker filtering is activated on the Oqus camera.
- There is a red warning sign  in the top left corner if the camera has problems calculating the markers. Hover the mouse over the warning sign to get more information.
The most common problem is that the camera does not have time to calculate all of the markers. This is actually caused by too much background light. If this happens, reduce the **Exposure time** or increase the **Marker threshold**.
- If the markers are grey and it says **Not used for tracking** in the middle of the view, that camera has been deactivated on the **Linearization** page in the **Project options** dialog.
- The delayed exposure setting is displayed next to the camera number, e.g. **expgrp: 1** when the camera is in exposure group 1.



3D view window and features

In 3D view windows the motion capture data is shown in 3D using the coordinate system that was used when calibrating.

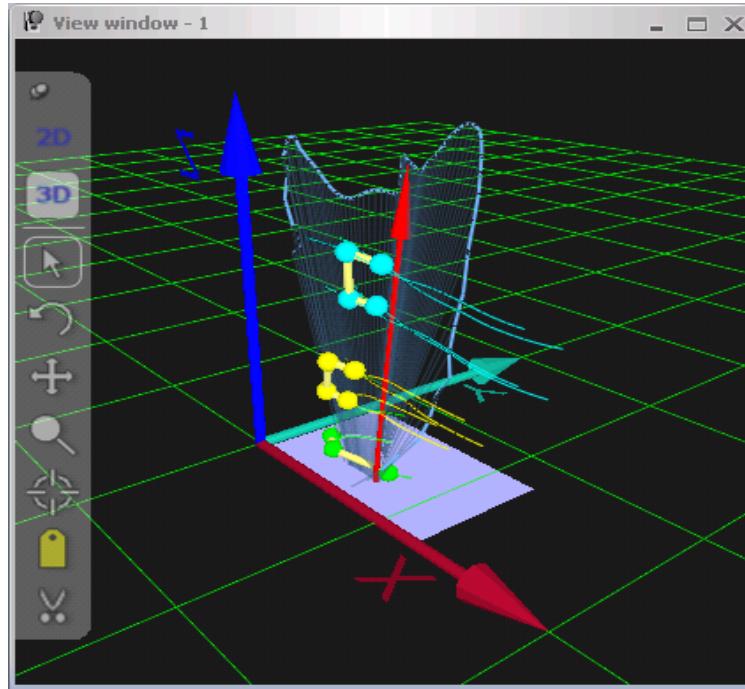

The 3D view window includes the following controls:

3D view toolbar

Contains tools for the 3D view, for example **Center trajectory** and **Quick identification**.

3D view window menu

Contains options for the 3D view and is opened by right-clicking in the 3D view window.

The following graphics can be shown in a 3D view window:

- Axes of the coordinate system of the motion capture (X = red, Y = light green and Z = dark blue).
- A green grid, which is the floor of the measurement volume (e.g. Z = 0, if Z is the vertical axis of the measurement setup).
- The position and orientation of each camera in this specific setup.
- Text labels for some of the graphics.
- Markers for the current position of the trajectories.
- Traces of the trajectories.
💡 Note: Only available in file mode.
- Bones between markers.
- 6DOF bodies.
- Purple squares, which display the location of the force plates.
- Red force vectors displaying the current force on the force plate.
- A blue force trace displaying the force applied to the force plate during the measurement.
💡 Note: Only available in file mode.
- Camera view cones.

- Covered and calibrated volumes.
- A white outlined box for the bounding box used by the 3D tracker.

All of these graphics can be toggled on and off and changed on the **3D view settings** page in **Project options** dialog.

It is also possible to change the viewpoint of the 3D view with the mouse. When actions are performed which change the appearance of the 3D view a red crosshair  is visible in the 3D view. The crosshair is always placed in the middle of the 3Dview window.

The following actions can be performed in 3D viewwindows to change the appearance of the 3D view:

Rotate the 3D view

Hold down the left mouse button in the 3Dview window and move the mouse, the view is rotated freely around the red crosshair.

Translate the 3D view


Hold down the right mouse button in the 3D viewwindow and move the mouse, the focus is moved in the 3D view.

Another way to move the focus is to hold down Shift+C and click on a trajectory then the focus is moved to that marker.

Zoom in/out the 3D view

Hold down both mouse buttons in the 3Dview window and move the mouse backward and forward. The view is zoomed out and in respectively, in reference to the position of the red crosshair. The mouse wheel can also be used to zoom. You can also zoom to the current position of the mouse cursor by holding down Ctrl.

Trajectories in 3D views

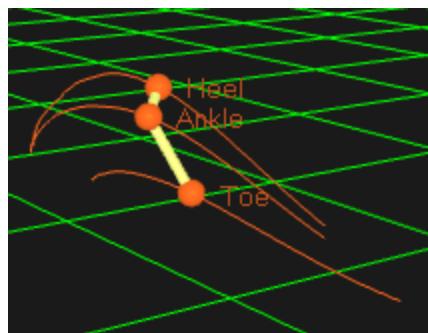

The trajectories of a measurement are shown with markers on their current position in the 3D view, where the colors of the trajectories are the same as in the **Trajectory info** windows. The data of the marker is shown in a **Trajectory info** window, see section “**Trajectory info windows and menu**” on page 28.

In file mode, the position throughout the measurement can be displayed with a trace, where the length of the traces is controlled by the bottom sliders on the **Timeline control** bar. You can change settings for the marker and trace display on the **3D view settings** page in the **Project options** dialog. Among other things you can turn on the display of the marker labels.

The following section describes some of the actions that can be performed on trajectories. For more information about the actions see the **QTM User manual**.

Drag and drop

The trajectories and traces can be drag and dropped to change their identity. For example drag and drop a trajectory from the 3D view to a label in the

Labeled trajectories window to identify it. If the label is not empty it is important that there is no part of the two trajectories that do not overlap completely. You can also drag and drop a trace on another trace to join them.

Quick identification

The quick identification method is the best way to manually identify the markers. First select an empty label. Then hold down Ctrl+Alt or select the **Quick identification** cursor to activate this method. Then click on a marker and it will be identified as the currently selected empty label in the **Labeled trajectories** window. For more information about quick identification see section “Manual identification of trajectories” on page 49.

Create bones

To create a bone, hold down Shift and select a pair of labeled trajectories by clicking in the 3D view window. Then press 'b' and there will be bone between the pair. Several bones can be created at the same time if the trajectories are selected one at a time with Shift + Left-click in the 3D view window.


Center on trajectory

Select a trajectory and press 'c' to change the viewpoint so that it centers on that marker. You can also use the **Center trajectory** cursor and click on the marker.




6DOF bodies in 3D views

Measured 6DOF bodies are displayed as separate local coordinate systems in a **3D view** window, where the origin is placed at the position of the local origin of the 6DOF body. The axes of the local coordinate system have the same color codes as the coordinate system of the motion capture.

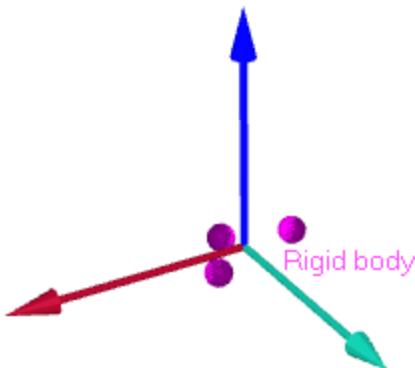

Virtual markers corresponding to the points of the 6DOF body definition are also displayed for each 6DOF body. The markers will therefore be displayed with two colors in the 3D view. The colors of the rigid body are set for each 6DOF body definition on the **6DOF tracking** page in the **Project option** dialog. The actual trajectories will automatically get a slightly brighter color.

Place the mouse over the local coordinate system to see information about the body in the current frame.

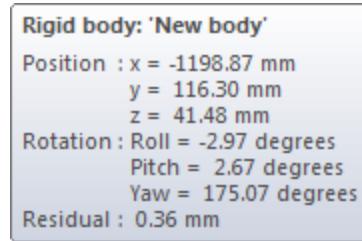

Trajectory info windows and menu

The trajectories of the measurement are listed in the **Trajectory info** windows. The following three **Trajectory info** windows are used in QTM:

Labeled trajectories window

Contains the trajectories that have been identified (labeled).

Unidentified trajectories window

Contains unidentified trajectories.

Note: After a motion capture the trajectories are placed by default in the **Unidentified trajectories** window, as long as identification is not performed in the automatic processing.

Discarded trajectories window

Contains trajectories that have been manually deleted.

By default the **Trajectory info** windows are placed on the right side of the main window. However they can be moved to any position.

The options for the trajectories are found on the **Trajectory info window** menu, which is opened by right-clicking on a trajectory or a group of selected trajectories. When multiple trajectories are selected the options on the menu is applied to all of the trajectories.



Data in Trajectory info windows

Labeled trajectories (18)							
Trajectory	Fill Level	Range	Type	X	Y	Z	Residual
+ r_shoulder	74.0%	183 - 700	Mixed	-652.81	-91.97	1483.13	2.64
+ th12	81.3%	95 - 663	Measured	-766.57	72.25	1125.80	1.73
+ r_axis	83.6%	116 - 700	Mixed	-586.68	-58.00	992.70	2.15
Part 1	76.6%	116 - 651	Measured	-586.68	-58.00	992.70	2.15
Part 2	0.4%	652 - 654	Gap-filled	-	-	-	-
Part 3	6.6%	655 - 700	Measured	-	-	-	-
+ r_suppat	81.7%	129 - 700	Measured	-429.41	-9.17	564.63	1.65
+ r_knjntln	85.6%	102 - 700	Mixed	-503.50	-69.44	488.78	1.14
+ r_tubtib	81.3%	132 - 700	Mixed	-430.83	-14.23	442.96	1.49
+ r_ankle	81.0%	134 - 700	Measured	-489.28	-30.62	55.64	0.54

The following data is listed in all of the **Trajectory info** windows:

Trajectory

The label of the trajectory and the color of its marker in 3D views. In addition to the color the symbol next to the label also shows if the trajectory and its trace are displayed in 3D views. The following symbols are used:

- Both trajectory and trace are displayed.
- Just the trajectory is displayed.
- , ● Nothing is displayed.

If a trajectory consists of more than one part there is plus sign (+) in front of the label. The parts are shown by clicking the plus sign or by using the left and right arrows.

 Note: In the **Labeled trajectories** window the label can be edited by double-clicking on the label. However, in the other two windows the labels are always **Unidentified** and **Discarded** followed by a sequential number, which can be used to separate the trajectories.

Type

The type of trajectory, which among others can be one of the following types.

Measured

A trajectory or a part that has been tracked from the measurement data.

Mixed

A trajectory with a mix of the other trajectory types.

Gap-filled

A trajectory or part that has been calculated with the gap fill function.

Fill level

The percentage of the current measurement range that the trajectory or part is visible. The measurement range is selected using the **Timeline control** bar.

Range

The range of frames, within the measurement range, with data for the trajectory or part.

X, Y and Z

The position (in mm) of the trajectory in the current frame. The coordinates use the coordinate system of the motion capture system set up when the system was calibrated.

Residual

The average of the different residuals of the 3D point. This is a quality check of the point's measured position.

 Note: The residual is not available for 2D tracked trajectories.

 Note: The number after the window name is the number of trajectories in the window.

Select and move trajectories

Select a trajectory in the **Trajectory info** window by clicking on the trajectory. It is possible to select multiple trajectories with standard Windows multiple selection methods. To deselect all trajectories, either click on the **Trajectory** heading to deselect all or use Ctrl+Left-click on the selected trajectories.

The trajectories can be moved to another **Trajectory info** window by a simple drag and drop, either one by one or in a group. If the trajectory is dropped in the **Labeled trajectories** window the trajectory is labeled 'New 0000', where the number depends on how many trajectories there are in the window.

Trajectories can also be joined by drag and drop. Select a trajectory and drag and drop it on another trajectory or on an empty label. The trajectory will then contain a new part. If the trajectories overlap there will be a clash error, which must be handled. Trajectories where one of the trajectories or a part is completely overlapped by the data of the other trajectory cannot be joined at all.

 Note: It is also possible to drag and drop trajectories from and within the 3D view window.

Data info window and menu

2d Data, Camera 2				
x	y	xSize	ySize	...
56644	36473	392	384	
66137	43471	488	512	
69986	46274	472	512	
58627	51787	490	576	

The data for the current frame can be viewed in the **Data info** window. Open this window by clicking **Data info window [Ctrl+D]** on the **View** menu. It is always placed on the left side of the main window. The following data can be viewed in the **Data info** window:

2D data for each camera


(Use the option **Current camera** if you want the 2D data of the last camera that you clicked on in a 2D view)

6DOF data with position and rotation

Analog data from the analog boards

Force data for the force plates

Use the **Data info window** menu to switch between the data types and to plot data, it is accessed by right-clicking in the **Data info** window.

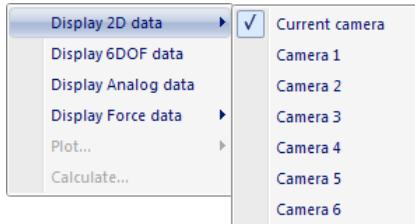


Plot windows


The **Plot** windows display a graph of the data that has been selected with the **Plot** command. The default window shows the graph and the **Legend** box, see screen dump below. The **Legend** box displays the names and colors of the series.

The graph contains the following default graphics:

- A title above the plot, which shows what is plotted in the graph.
- A purple line in the plot, which represents the current frame in the **View** windows.
- An X-axis, which is always in frames.
- A Y-axis, where the unit depends on what is plotted in the window.

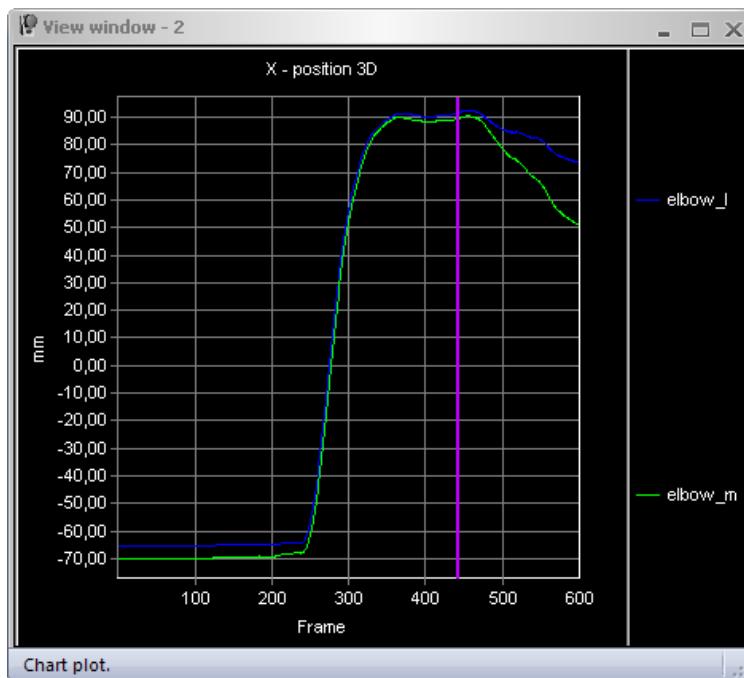

The settings for the **Plot** window are reached by right-clicking in the plot area, which will open the following menu.

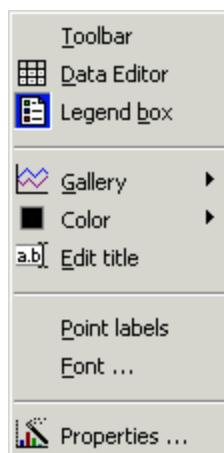

Open the toolbar to get more settings, including a zoom tool.


Project options dialog

The configuration of software and hardware settings in QTM is accessed via the **Project options** on the **Tools** menu. The settings are described in detail in the **QTM User manual** and some basic settings for a capture are described in the chapter “QTM capture settings” on page 33.

 Note: Most of the settings will only affect a measurement if QTM is still in the preview mode, i.e. before the measurement has started.


QTM capture settings

Capture settings overview

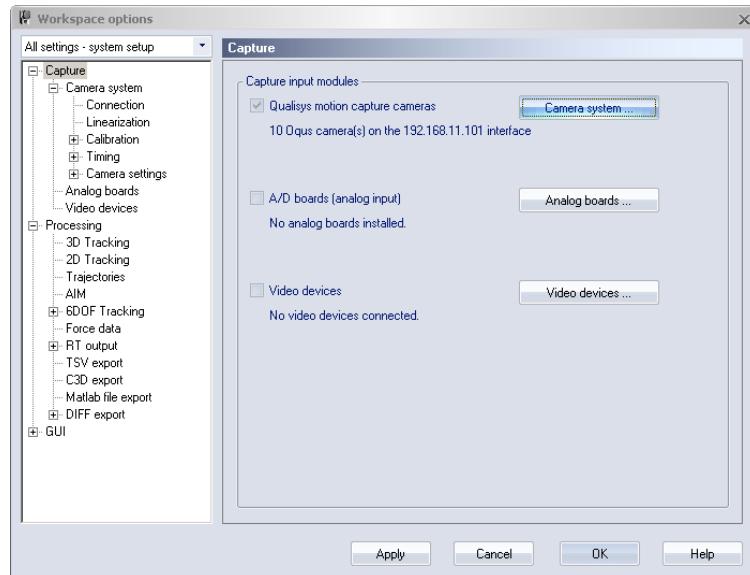

The **Capture** tree in the **Project options** dialog contains all of the settings used by the camera system and hardware in a measurement. Only the most important settings are described below, for a detailed description please refer to the **Qualisys User manual**.

Connecting to the system

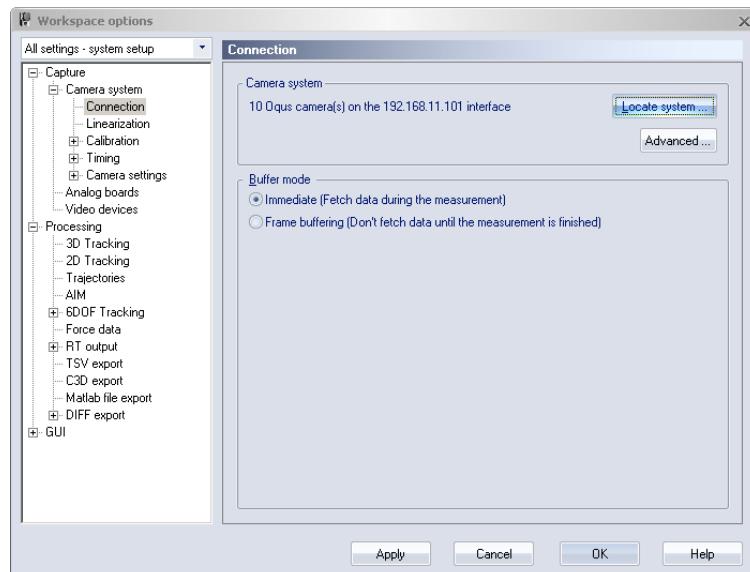


Before the first use of the camera system it must be connected to QTM. Make sure that the camera system hardware has been correctly installed.

The connection of the camera system to the QTM software is performed on the **Connection** page in the **Project options** dialog. When the connection has been made, the camera system can be used as long as the properties are not changed.

Follow these steps to connect the camera system to QTM:

1. Switch on the camera system, wait for the cameras to start up properly and start QTM.
2. Open the **Project options** dialog and go to the **Connection** page.
3. Click **Locate system**.
4. Click **Start**.
5. Choose the camera system and click **OK**.

Loading linearization files

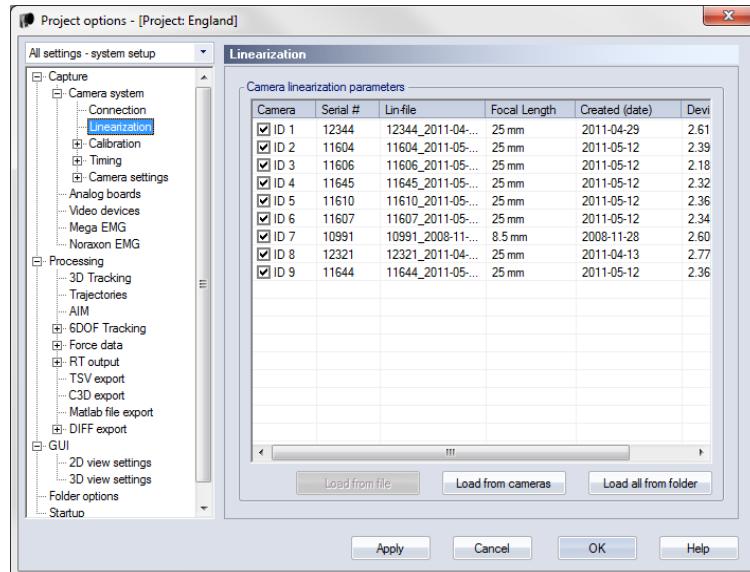


A camera transmits coordinates in two dimensions for each marker's center point in the field of view of the camera. The locations of the center points reflect the optical view of the sensor chip in the camera. Since this cannot be done without distortion, the QTM software must have information about the distortion of each camera in order to get correct 3D data. This is achieved by an adjustment of the camera data with a linearization procedure in QTM, which is then saved in a linearization file and stored in the camera.

The linearization files of the current system are specified on the **Linearization** page in the **Project options** dialog. To load the linearization files stored in the cameras click on **Load from cameras**. A dialog appears and QTM starts to load the linearization files from the Oqus cameras.

Calibration settings

To receive 3D data the system must be calibrated. The most commonly used calibration is the Wand calibration that uses two calibration objects to calibrate the system: the L-shaped reference structure and the calibration wand.

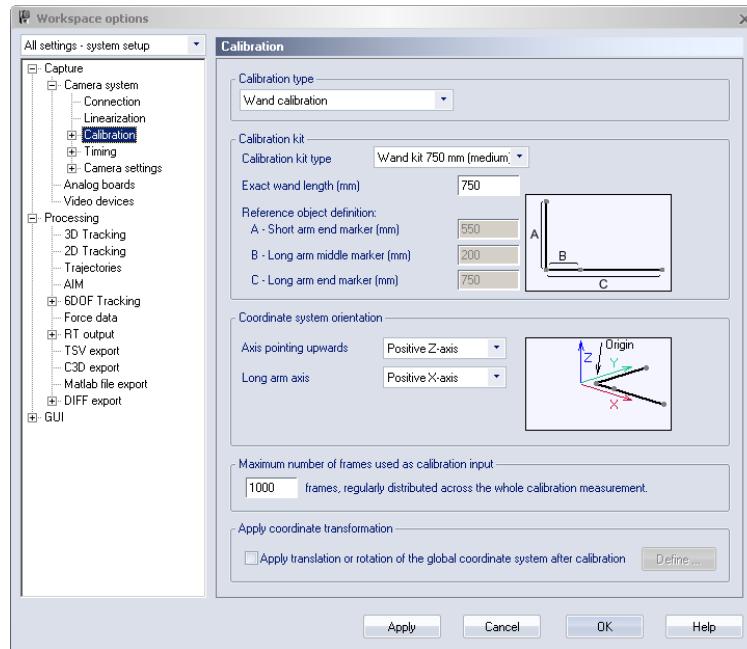

The settings that must be specified for the calibration are the **Calibration kit type** and the **Exact wand length**. Without these settings the calibration will be incorrect.



Oqus marker settings

The two important settings for marker calculation are **Exposure time** and **Marker threshold**. Usually it is sufficient to change them to see the markers in QTM. However it is not possible to give exact advice on how to set the settings, because the relations are too complex. The tips below will show you the basic method of how to use the settings. There are also other settings that can help to get a better measurement, e.g. **Marker masking** and **Marker discrimination**.

First of all make sure that the focus and aperture are correct, see section “Setting the aperture and focus” on page 19.

Exposure and threshold must be used together, because changing one can lead to that you have to change the other. Follow these steps to set exposure and threshold settings.

1. Start a new measurement and stay in RT/preview mode. Use the **Camera settings** sidebar to change the settings.



2. Set the **Exposure time** to 400 microseconds (200 μ s for Oqus 3+ and Oqus 4) and **Marker Threshold** to 175. These are the default values and often it is a good starting point.
3. Start with looking at the marker intensity image. Click on the **Marker intensity** button in the **Camera settings** sidebar to show the marker intensity image.
4. If the markers are not bright red in the marker intensity image, try increasing the **Exposure time** until they are bright red. For example with a longer distance to the markers, you may need a longer exposure time.
 - a. For higher frequencies it might not be possible to increase the exposure so that the markers are bright red. However, as long as they are brighter than the background the camera should be able to find them by lowering the marker threshold.
 - b. If there are extra reflections, you can try reducing the **Exposure time**. Extra reflections are anything that is not a marker and has a color different from blue. Green is the threshold level.
 - c. Remember that you can set this setting individually for each camera.
5. Switch back to Marker mode and check if the markers are visible. If they are not, go back to Marker intensity mode and adjust the **Marker threshold** value. It is not possible to give an exact value for the threshold at a certain exposure, because each setup is different.
 - a. Increase the threshold if the background looks light blue or even green in the marker intensity image. A light blue background will make it harder for the camera to calculate the markers.
 - b. Decrease the threshold if the markers are not bright red in the marker intensity image. For example, at short exposure times of 200 μ s (100 μ s for Oqus 3+ and Oqus 4) and lower the threshold needs to be low, usually around 100-150.
 - c. Make sure that the marker calculation is stable at the selected threshold. If the threshold gets too low there will be a lot of extra reflections or no markers at all. A too high threshold will result in small markers or missing markers.
 - d. Remember that you can set this setting individually for each camera.
6. Finally check that the markers are large enough to give a good accuracy. Check that the marker size is at least 200 in the **Data info** window.

Calibration

Calibration overview

The QTM software must have information about the orientation and position of each camera in order to track and perform calculations on the 2D data into 3D data. The calibration is done by a well-defined measurement procedure in QTM.

There are three methods that can be used to calibrate a camera system: Wand calibration, Frame calibration and Fixed camera calibration. Only Wand calibration is described in this Getting started manual.

The following items are important to think about before the calibration:

- Before the calibration, make sure that the calibration settings are correct on the **Calibration** page in the **Project options** dialog. The focus and aperture of the cameras must also be set before the calibration.
- When using Wand calibration or Frame calibration it is important that the camera system has been placed correctly to achieve a high-quality calibration, see section “Calibration overview”. For Fixed camera calibration it is equally important with the camera positions but in this case the cameras are fixed and must be installed correctly before the calibration.
- It is recommended that the motion capture system is calibrated before each measurement session to make sure that the captured data has high quality. Check regularly that the calibration is OK during long measurement sessions.

There is an automatic check if the cameras has moved, which is activated on the **Calibration Quality** page in the **Project options** dialog.

 **Important:** Each time a camera in the system is moved (even the slightest) a calibration must be performed before high-quality data can be captured again. However, some errors with a calibration can be fixed with reprocessing.

Camera positioning

Cameras must be mounted firmly on tripods or other stable structures, which isolate the camera from movements or vibrations of any sort. Once the system has been properly set up, the cameras must be arranged to fit the measurement setup. When arranging the cameras, it is best if they are in operating mode.

1. Start the measurement computer and the QTM software.
 2. Open a new file with a **2D view** window.
 3. Arrange the cameras to cover the entire measurement volume of the wanted motion. It is suggested to mark the corners of the measurement volume with markers, and use the **2D view** window to make sure that each camera can see the markers as expected. The **Video** view can also be used to see the camera field of view.
 - The best setup for a 3D motion capture system is to position it so that all cameras can see the L-shaped reference structure during the calibration.
-  Note: The cameras can be positioned so that just two of the cameras are able to see the calibration reference object. The rest of

the cameras must then overlap each other's field of view (FOV) to be able to calibrate the system. For this setup QTM will automatically use the Extended calibration method.

- To reconstruct 3D data, at least two cameras must see each marker during the measurement. Therefore it is best to position the cameras so that as many cameras as possible see each marker during the measurement.
- The angle of incidence between any two cameras should ideally be more than 60 degrees and at least more than 30 degrees. The accuracy of the 3D data calculated from only two cameras placed at less than 30 degrees can degrade below usable levels.
- In order to avoid unwanted reflections, position the cameras so that every camera's view of flashes from other cameras is minimized. E.g. put the cameras above the measurement volume so that the cameras have an angle of about 20 degrees in relation to the floor.


How to perform a Wand calibration

The Wand calibration method uses a calibration kit that consists of two parts: an L-shaped reference structure and a calibration wand.

Place the L-shaped reference structure so that the desired coordinate system of the motion capture is obtained. It is best if all cameras in the system can see all markers on the reference structure. This is not an unconditional requirement but it gives the highest accuracy. If some cameras cannot see the reference structure the Extended calibration method is used automatically.

The calibration wand is moved inside the measurement volume in all three directions. This is to assure that all axes are properly scaled. The calibration algorithms will extract each camera's position and orientation by evaluating the camera's view of the wand during the calibration. To achieve a good calibration of the system, it is recommended that you follow these instructions on how to move the wand.

- To assure that all cameras can see the wand as much as possible during the calibration, make sure that the cameras are not blocked by the person moving the wand around in the measurement volume. This ensures that no camera will be blocked for a longer time period.
- One suggestion on how to move the wand, is to move it in one direction at a time. Start by holding the wand positioned in the Z direction, i.e. the straight line between the two wand markers should be parallel to the Z axis. Move the wand in the entire measurement volume. It is important to fill the entire measurement volume with calibration points. Make sure that both the lower and upper parts of the volume are covered. Repeat the same procedure with the wand positioned in the X and Y direction. It is particularly important to collect points where there will be many markers during the motion capture.

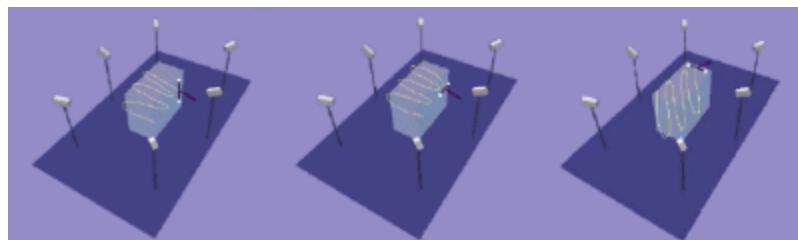

 Note: In this picture the reference structure is not indicated to make the picture more distinct. The reference structure must of course always be present during the calibration. The box in the figure represents the measurement volume.

There are many other ways to move the wand and the optimal moving method varies depending on the application, i.e. you have to find out the best method for your application. It is recommended that the moving method is systematic and easy to repeat.

The steps below are an outline of what should be done to calibrate the camera system with Wand calibration method.

Follow these steps to calibrate the camera system:

1. Switch on the camera system and start QTM.
2. Open a new file by clicking the **New file** icon .
3. Place the L-shaped reference structure in the measurement volume.
4. Set the settings on the **Calibration** page in the **Project options** dialog.
5. Click **OK**.
6. Click the **Calibration** icon  and set the settings in the **Calibration** dialog.
 Note: A suitable setup to assure a correct calibration of all cameras is to use 20-30 seconds for the calibration. It is, however, very dependent on the size of the measurement volume. You need to specify enough time to cover the entire volume with the wand.
7. Click **OK**.
8. Move the calibration wand in the measurement volume.
9. Check that the calibration is OK in the **Calibration results** dialog and click **OK**.

 Note: Each calibration is saved in a QCA file in the **Calibration** folder in the **Qualisys Track Manager** folder. The file name contains the date and time of the calibration.

Capture

Capture overview

QTM can capture different data including marker, video, force, EMG and analog data. There are two ways of collecting data in QTM, either capture a file in QTM or run QTM in RT (preview) mode and then collect the data in another program. Which one you use depends on your application.



Real time mode

The QTM real time process enables QTM to send data to any program that can receive data on TCP/IP or UDP/IP. At the moment there are interfaces for iOS devices, Visual3D, Motion Monitor, Motionbuilder, LabView and Matlab. Please contact Qualisys AB if you want to use these interfaces. If you want to transfer real time data into other programs please contact Qualisys AB. You will then be provided with the real time protocol, so that you can write your own interface for that particular program.

QTM will process the 2D in real time as soon as a measurement file is open in preview, i.e. when you click on **New** in the **File** menu. The data can then either be displayed in QTM or sent to another program.

The real time performance depends on the computer, the number of cameras and the number of markers. To get the best performance, use a dual-core desktop computer with at least 2 gigabyte RAM and a good graphics card.

The following processing steps can be done in real-time, depending on the settings on the **Processing** page in the **Project options** dialog.

3D tracking

Apply the current AIM models

Calculate 6DOF

Calculate force data

This data can then be viewed in QTM and can also be sent to another computer via TCP/IP. These are the different steps in the real-time process.

1. The cameras capture data at a fixed frequency.
2. RT processing of data in QTM which is performed as fast as possible.
3. Data is sent on the TCP/IP output. If QTM takes too long time processing a frame, the next frame captured by the cameras will not be processed at all. This is shown in the main status bar as decreased a RT frequency.


How to use real time

The following steps describe how to use real-time and send the data to another program, for example Visual3D.


1. Before you start the real-time, make sure that you have an AIM model or 6DOF bodies for the movement that you are going to capture.

 Note: It is best if the AIM model has been created on a measurement made on the subject that will be captured in real time.

2. Switch on the camera system and start QTM.
3. Open a new file by clicking the **New file** icon .
4. Open **Project options** dialog and go to the **Processing** page.
5. Activate the **Real-time actions** that you want to use. For example for Visual3D the following actions should probably be used: **Track each frame: 3D**, **Apply the current AIM models** and **Calculate force data**.

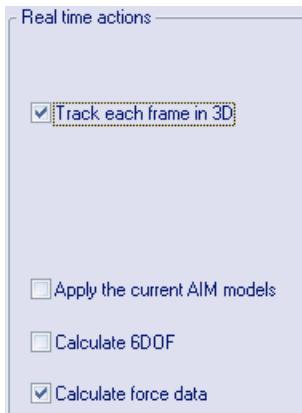

6. Check the settings for the actions that have been selected.
7. Go to the **Camera system** page and set the capture rate. You can either use the same frequency as the **Marker capture frequency** or set a lower frequency with the **Reduced real time frequency** option. The maximum capture rate depends a lot on the computer and the AIM model, but 100 Hz usually works with a dual-core computer.

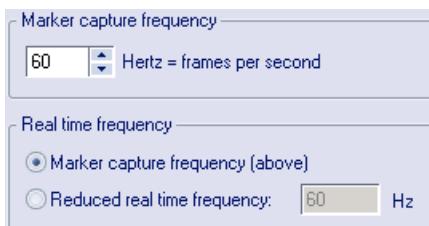

8. Go to the **GUI** page and set the **Real time mode** screen update to 15 Hz. This is only to reduce the processing needed for the display update.
9. Test the real-time with the motion that you want to capture. Look especially at how the AIM model is applied and if the RT frequency shown in the **Status** bar is close to the capture rate. If it differs too much, lower the capture rate, it might also help with a new AIM model or changing the tracking parameters.
 - a. If the real-time is slow close all windows, including the **Data info** window and the **Trajectory info** windows, except for a **3D view** window. Activate only the necessary processing steps on the **Processing** page.
 - b. When the real-time is working fine you can even turn off the GUI with the **Disable GUI update** button . This will reduce the processing capacity needed by QTM.
 - c. AIM can be restarted with the **Apply AIM model** button .
10. When you are satisfied with the real-time in QTM, you can connect to QTM from the other program.

Capture measurement files


The reason for capturing data is to store the data from the cameras. Captured data can then be processed and exported to other programs. The QTM files store all of the data that is captured, including marker, analog and force data. Video files can also be captured with the Oqus camera or a DV camera, but are stored as separate AVI files.

How to capture a file

Starting a capture is simple, follow the steps below.

1. Switch on the camera system and start QTM.
2. Before you can start a measurement you must open a new empty capture file with **New** on the **File** menu or by clicking the **New file** icon . The file will be opened in preview mode where you can check measurement volume and settings before starting the capture.
3. Check that the system is calibrated.
4. Check that the markers are visible, otherwise change the camera settings see chapter “Oqus marker settings” on page 35.
5. It is a good idea to check the settings on the **Processing** page in the **Project options** dialog, see chapter “Processing overview” on page 47. There you can define which processing steps will be automatically applied to the file.
6. Click the **Capture** icon .
7. Specify the capture settings in the **Start capture** dialog.

Capture

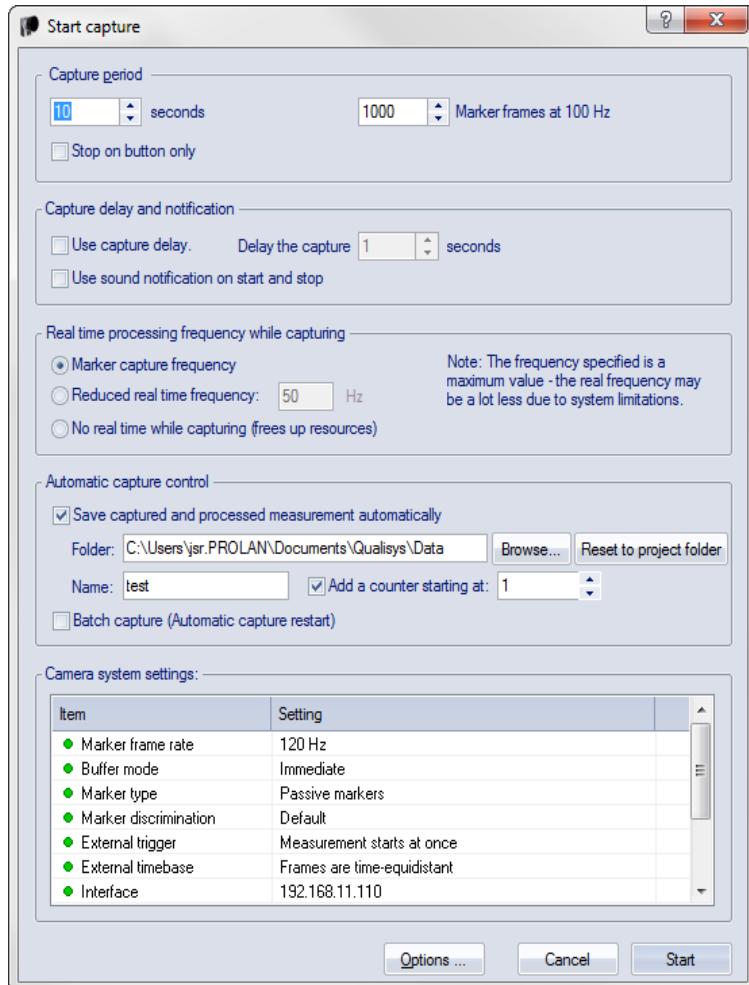

Check especially the following settings:

Capture period is how long the capture will be. You can also see the capture rate and the number of frames to the right of the setting.

Marker real time frequency while capturing controls the preview/RT processing during the capture. When this is turned on, RT data can be retrieved by another program during the measurement. If it is turned off, there will be no preview in QTM and no RT data can be retrieved by another program.

The **Automatic capture control** settings can be used for saving the files automatically and also for activating batch capture.

8. Click **Start** to start the capture. When the capture is finished, all of the processing steps will be performed and then the new motion capture file is displayed in QTM, unless batch capture is activated. In batch capture, QTM will immediately start waiting for the next capture as soon as the processing of the previous capture is finished.

The capture can always be stopped with **Stop capture** on the **Capture** menu or by clicking the **Stop capture** button ■■■■■.



How to capture an Oqus video file

The following outline describes how to capture high-speed video with an Oqus camera.

1. Switch on the camera system and start QTM.
2. Open a new file by clicking the **New file** icon  If you want to use the 3D overlay functionality it is important to check that the camera system has been calibrated.
3. Switch to 2D view if it is not visible, i.e. right-click in the View window and select **Switch to 2D view** from the menu.
4. Right-click on the 2D view for the camera that will capture high-speed video. Select **Mode/Video** to switch to Video mode. The 2D view for that camera will switch to a video image.
5. Open the aperture and set the focus.
6. Change the settings for Video capture in the **Camera settings** sidebar.
 - a. Set the video **Capture rate**. This can be set independently of the marker image rate.
 Note: The image size is reduced automatically when it is changed on the **Camera settings** sidebar.
 - b. Set the **Exposure time** to a value that makes the image bright enough, test until you are satisfied.
 If you have no extra light the exposure time needs to be quite high, at least 16000 microseconds or even up to 40000. This limits the capture rate that can be used. It also means that fast movements will be blurred.
 For high capture rates and measurements with fast movement, extra lighting is needed because the exposure time must be extremely short, sometimes as short as 100 microseconds.
7. Use a Codec if you want to reduce the size of the avi-files. If the file will be used in analysis it must not be reduced too much as it can influence the analysis.
 Note: The codec does not change the download speed.
 Note: The **Image format** will reduce the image directly in the camera, but then you will lose pixels.
8. A video capture takes about a second to initialize so it is recommended that you use external trigger and pretrigger to start the capture.
9. Click the **Capture** icon  and set the capture time and other capture settings. Keep the capture time as short as possible to reduce the fetching time.
10. Click OK and then press the trigger button when you want to start the measurement.
11. Wait for the fetching of video data to finish. Depending on how large file you have captured and the number of cameras, it can take up to some minutes to fetch the data.
12. The video file is displayed in the 2D view window and can be played in synchronization with other data.

 Note: If the video data is uncompressed it is not possible to play the file in normal speed as the video playback will be too slow. Then it is better to step through the file instead.

Batch capture

With batch capture, QTM will capture several consecutive measurements. Batch capture is activated with the **Batch capture** option on the **Start capture** dialog, before the start of the measurements. In this dialog, the options **Save captured and processed measurement automatically** and the automatic counter must also be selected so that each measurement is saved in a separate file.


Before each new measurement in a batch capture QTM will wait for a start signal and the whole batch capture is stopped by a stop signal. These signals are given in different ways depending on whether external trigger is used or not.

External trigger

If the external trigger is used to start each measurement, to stop the batch capture you need to press Esc or click **Close** on the **File** menu. **Stop capture** on the **Capture** menu can be used during an individual capture to stop just that capture.

No external trigger

Start each measurement by clicking **Yes** in the **Next measurement** dialog. Stop the batch capture by clicking **No** in the dialog. **Stop capture** on the **Capture** menu can be used during an individual capture to stop just that capture.

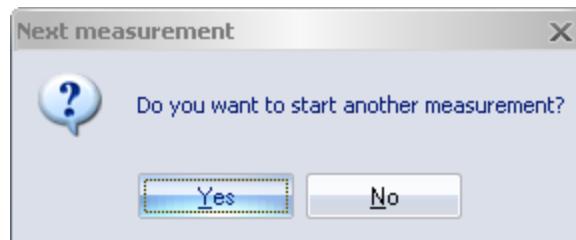

 Note: All of the processing steps that are selected on the **Processing** page in the **Project options** dialog will be performed before the next measurement can start.

Processing

Processing overview

After the capture is finished the data must be processed in QTM to get 3D data, 6DOF data and force data. Some functions in QTM, such as AIM and exports, can also be performed in processing. The processing can be done automatically after the capture, manually in the current capture file or automatically on saved files by using **Batch processing**.

The settings of the automatic processing, which are the same for the processing directly after capture and batch processing, are specified on the **Processing** page in the **Project options** dialog. The following processing steps are available on the page:

1. **Auto backup**
This option will create an automatic backup file containing the captured 2D data before starting the other processing steps.
2. **Track the measurement**
 - a. **3D**
 - b. **2D**
3. **Merge with Twin slave**
4. **Gap-fill the gaps**
5. **Apply the current AIM model**
6. **Calculate 6DOF**
7. **Calculate force data**
8. **Export to TSV**
9. **Export to C3D format**
10. **Export directly into Matlab**
11. **Export to Matlab file**
12. **Export to DIFF format**

The settings of each processing step can be reached from the **Processing** page, except for **Export directly into Matlab** which has no settings. This manual only describes settings for AIM. To read more about the settings for other steps, please refer to the **QTM User manual**.


Batch processing

With batch processing, several capture files can be processed with the same settings. The same actions can be performed in batch processing as in the processing that is performed after a measurement has been captured.

 **Important:** If **Track the measurement** is used when batch processing, the tracking of the saved files will be changed. This means that all of the trajectories will be unidentified again unless you apply AIM in the same batch process.

To start a batch process, follow the steps below.

1. Click **Batch process** on the **File** menu to open the **Batch processing** dialog.

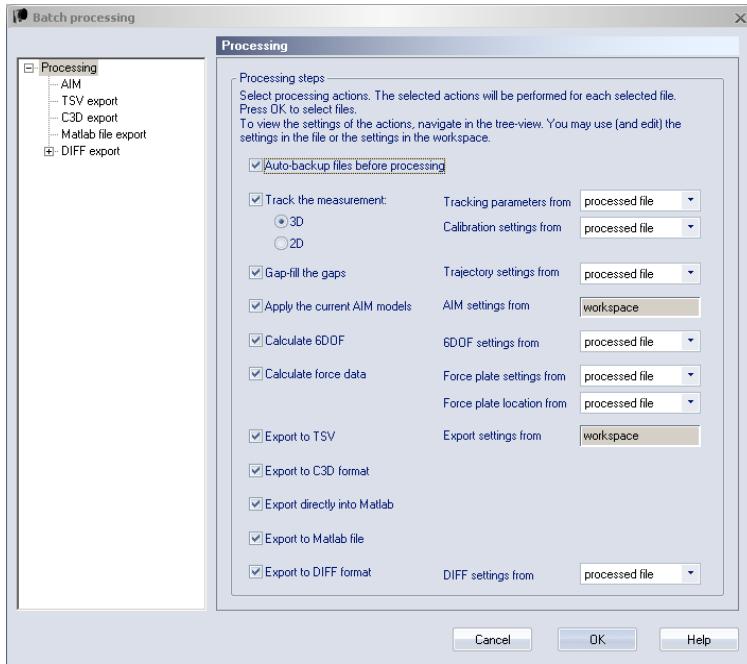

2. Select the processing steps that you want to perform from the list. The processing steps are the same as those that can be applied in post-processing.
3. Choose the source for the settings from the drop-down lists to the right. For some settings the only possible source is the workspace, for the other settings there can be three possible options:

Processed file

The settings used are the ones present in each processed file and cannot be edited. This option is always selected by default because the other options remove the original settings from each file.

Project

The settings are copied from the current workspace to each processed file. The settings can be edited in the tree-view to the left. Editing the settings will change the current workspace settings as well. This option is often the best to use when you want to change the settings of the processed files.

Present file

The settings are copied from the present file (the file open in the QTM window where you opened the **Batch processing** dialog) and can be edited in the tree-view to the left. Only available for **Calculate 6DOF**, **Calculate force data** and **Export to DIFF format**.

Note: For **Track the measurement** the settings are split into **Tracking parameters** and **Calibration settings** so that the sources of these settings can be set separately. Make sure that the calibration is valid for all of the files you want to batch process when using the option **workspace**.

Note: For **Calculate force data** the settings are split into **Force plate settings** and **Force plate location** so that the sources of these settings can be set separately. This can for example be used on files where the

location varies, but you want to process the files with the same **Force plate settings**. In this case, select **processed file** as the source of **Force plate location** and select **workspace** as the source of the **Force plate settings**. Edit the force plate settings by selecting them in the tree on the left.

4. Click **OK**. Select the files and click **Open**. The actions are performed on the files and then the files are saved and closed.

 Note: Any changes to a capture file are saved automatically, i.e. it will be saved if any step except the export options are used.

Identification of trajectories

The first step when working with the data is to identify the trajectories by giving them the correct label. This is either done manually or automatically, with the use of an AIM model. When processing a lot of measurements with the same marker setup it is recommended to use an AIM model as it will significantly reduce the processing time.

Manual identification of trajectories

Even if you plan on using an AIM model to identify the trajectories, you need to identify at least one file by hand, since an identified file is needed to create an AIM model. Manual identification can be performed with any of the following three methods, and they can also be used together.

Quick identification method

Manual identification is best done with the quick identification method. Follow these steps to use this method:

1. To use the quick identification method you need to have a label list in the **Labeled trajectories** window. Load a label list or enter the labels manually with **Add new label** on the **Trajectory info window** menu.
2. Select the first label in the **Labeled trajectories** window that you want to identify.
3. Hold down Ctrl + Alt or select the **Quick identification** cursor  and click on the marker or trace in the **3D view** window that matches the selected label.
 - a. When you click on the marker or trace, it will be added to the selected label. If you have not selected any label, the marker you click on will be added to the first empty label in the list. If there is no empty label in the list, the marker you click on will be added as a new label at the end of the list, and you can edit its name at once.
 - b. You can also add additional parts to a trajectory with the quick identification method. Select the labeled trajectory in the list, then hold down Ctrl + Alt and click on the unidentified marker. However, if the two trajectories cannot be joined (because of too much overlap) the selection will just move to the next label without any action.
 - c. By holding down Shift when clicking on the trajectory you can join a trajectory to the previous label in the **Labeled trajectories** window. This option can for example be used if the most recently identified trajectory only covers a part of the measurement range. In

this case, you can add more parts by clicking on the next part while holding down Shift.

4. QTM will move the selection automatically to the next label in the list. Hold down Ctrl + Alt every time you click on the next corresponding marker. Continue until all markers in the **Labeled trajectories** window are identified.
5. You can start and restart the process on any label in the list.
 - a. This means that if you click on the wrong marker, you can undo the identification with Ctrl + z and then make sure that the correct label is selected and use quick identification again.

Drag and drop method

Manual identification can also be performed by dragging and dropping trajectories. The following drag and drop methods are available to identify trajectories.

From the **Unidentified trajectories** window to the **Labeled trajectories** window.

From the **3D view** window to the **Labeled trajectories** window.

One trajectory (trace) in the **3D view** window to another trajectory (trace) in the **3D view** window.

A trajectory can be dropped either on a label to add the data to that label or in the empty part of the **Labeled trajectories** window to create a new label. If it is dropped on a label that already contains data, the two parts will be joined if they do not overlap too much.

Identify method

It is also possible to use the **Identify** option on the **Trajectory info window** menu, which appears whenever you right-click a trajectory, either in a 3D view window or in a Trajectory info window.

There are two more features that are useful when identifying trajectories.

- The keyboard shortcut **c** can be used to center on the selected marker in a **3D view** window.
- The **Trajectory info window** menu option **Center trajectory in 3D**, will also center on the selected trajectory or part. However, if the trajectory is not visible at current frame it will also move the current frame to the first frame of the trajectory.



Automatic Identification of Markers (AIM)

Automatic identification of trajectories in the QTM software is performed by a module called AIM (Automatic Identification of Markers). The AIM model is created from identified files and can then be applied to any measurement that captures similar motions with the same marker set. This makes AIM very flexible and powerful - any motion can be made into an AIM model.

Generating an AIM model

AIM needs a model file to work, which is generated from a measurement or multiple measurements with identified trajectories.

⚠️ IMPORTANT: To get the best AIM model it is important that you continue to add QTM files with the same marker setup to the AIM model, if the automatic identification fails. You can add files with subjects of different sizes, as long as you have used the same method to place the markers. Use the **Add to existing model(s)** setting when you identified the file manually, because if you create a new model for each subject you will not use the full potential of the AIM functionality.


Follow the steps below to add a file to a model file or generate a new model file.

1. Make a measurement and identify all of the trajectories that will be included in the model. The model becomes better if you follow the guidelines in “Guidelines for data added to AIM models” on page 53: The two most important guidelines are:
 - Make sure that the subject in the file is moving. The movement does not have to be exactly the same as the rest of the captures, just that it includes most of the movement. When you add more measurements to an AIM model, it becomes less important that the added files include a lot of movement.


 Note: Even if you want to make an AIM model for a static measurement it is important to have some movement in the first file. This is to make sure that the internal bone definitions are correct. However, the subsequent files can be static because then the definition is already correct.

 - Make sure that the trajectories have the correct identity throughout the file. You can select a smaller measurement range to delete the incorrect parts if you do not want to identify the whole file.

2. Select the measurement range of frames that will be used in the model. Use the scroll boxes  on the **Timeline control** bar to select the range. Choose a range where the motion is typical for what you are going to capture, most of the time it is best to use the whole measurement range unless there are large gaps or incorrect data.

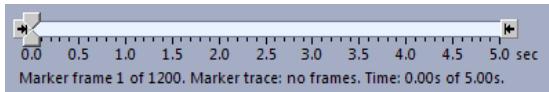

3. Click the **Generate model** icon  on the AIM toolbar or click **Generate model** on the **Auto ID** menu.
4. Select whether you want to **Create new model** or **Add to existing model**. Use add to existing model if you already have an AIM model with the current marker setup.

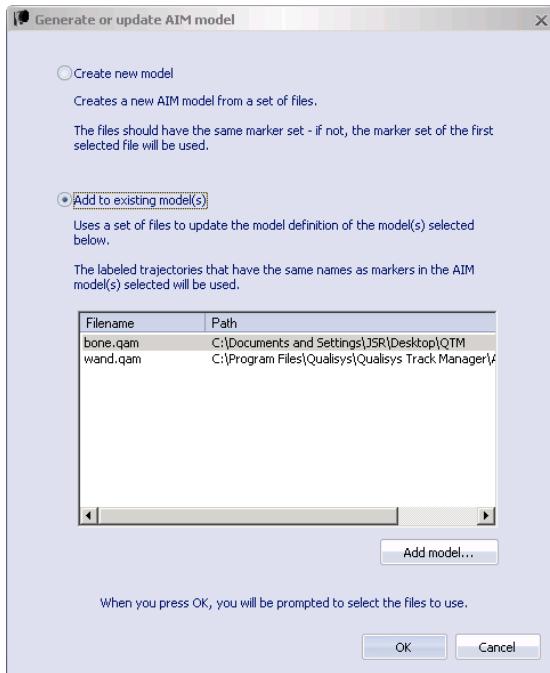

Create new model


Select the directory where you want to save the model, e.g. the same folder as the capture files. Then enter a **File name** for the model and click OK. The model is generated automatically.

Add to existing model

Select a model or models from the list to which you want to add the movement in the current file. This option is used to extend the motion of an existing model and therefore improve how it is applied to a new subject. It is the recommended option if you already have a model with the same marker setup. The previous AIM model is automatically saved in the AIM models/Backup folder.

- The models in the list are the models that are available on the **AIM** page. By default all the models in the **Applied models** list are selected. Click Add model to browse for another model.
- The file must include all the labeled trajectories that are included in the model. It can however have other labeled trajectories as well, only those with the same name as in the model will be used to update the model. This means that you can update several models at the same time. E.g. if the file include three person with different names on their respective labeled trajectories, then you can select all three AIM models in the list at once and all three models will be updated.

5. The dialog below is displayed when the AIM model has been created or updated. If anything goes wrong in the process it will be displayed in the **Results** list. For example you will get a warning if the first file is static or if the added file is too different and in those cases you can choose to abort the process.

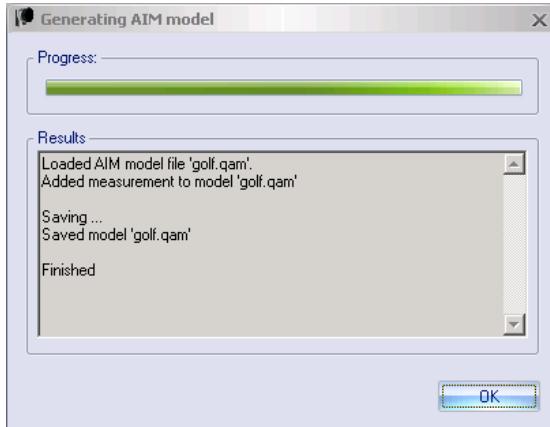

If you created a new AIM model it is then loaded to the **Applied models** list on the **AIM** page in the **Project options** dialog. Any other model in this list will be moved to **Previously used models** list.

If the data is added to an existing model, then the model files are updated and nothing is changed on the **Applied models** or **Previously used models** lists.

The model can now be applied to other measurements.

 Note: The procedure is the same if you use the **Generate from multiple files** option on the **AIM** menu, except that you select the files you want to use in a dialog. However it is important that the files have the same labeled trajectories.

Guidelines for data added to AIM models

Follow these guidelines when you check the data before adding it to or creating an AIM model.

- Make sure that the subject in the file is moving. The movement does not have to be exactly the same as the rest of the captures, just that it includes most of the movement. When you add more measurements to an AIM model, it becomes less important that the added files include a lot of movement.

 Note: Even if you want to make an AIM model for a static measurement it is important to have some movement in the first file. This is to make sure that the internal bone definitions are correct. However, the subsequent files can be static because then the definition is already correct.

- Make sure that the trajectories have the correct identity throughout the file. You can select a smaller measurement range to delete the incorrect parts if you do not want to identify the whole file.
- The colors of the trajectories and any bones between them are also saved in the model. For example to make it easier to verify the identification after AIM has been applied, the colors of the trajectories can be set with **Set different colors** on the **Trajectory info window** menu before creating the AIM model.
- Trajectories that are left unidentified or discarded will not be included in the model.
- If you have several subjects it is recommended to make an AIM model for each subject, see “AIM models for multiple subjects” in the QTM user manual.

The following two steps are not as important when you add files to the AIM model.

- The AIM model will be better if all of the trajectories are visible throughout the whole measurement. Therefore it is a good idea to gap-fill the trajectories as long as the gaps are relatively small. However if there are large gaps it is sometimes better to omit the gaps by limiting the measurement range in step 2 below. This is especially important if you use clusters, because then you have several trajectories that move along the same path.
- The data of each trajectory should be as good as possible. This is especially important if your model includes trajectories that are close to each other. Then a small erratic noise on a trajectory can make two trajectory paths come very close to each other, which makes the identification difficult. Therefore if the data is erratic, i.e. the trajectory movement is not smooth when you play the file, you should delete the erratic data by the following process:
 - Find the frame where the erratic data starts. Split the trajectory at that frame.
 - Then step through the frames to locate where the trajectory data is OK again. Split the trajectory again.
 - Delete the part with erratic data that you have just created.
 - Repeat these steps for all of the frames where you can find erratic data and then gap-fill the trajectories


Applying an AIM model

The AIM model can be applied to files with captured motions that are similar to any part of the motion in the model. I.e. if your model includes a lot of different motions made by a human, then the captured trajectories can be identified if another human makes one or more of these motions. An AIM model can be applied either as a processing step or manually.


Below follows a description of how to apply an AIM model manually on a capture file.

1. Open a capture file in QTM.
2. Make sure that all of the trajectories are in the **Unidentified trajectories** window or the **Identified trajectories** window. Discarded trajectories are not used by AIM. It is also important that all required parts of the measurement is included in the selected measurement range, since the AIM model is only applied to trajectories with parts within the selected measurement range.
3. Click the **Apply model** icon  on the AIM toolbar or click **Apply model** on the AIM menu. The **AIM application settings** dialog is displayed.


 Note: The AIM settings are copied from Project options, if you want to use the AIM settings and models saved in the file you need to reprocess the file and select from **measurement** as source, see chapter “Processing overview” on page 47.

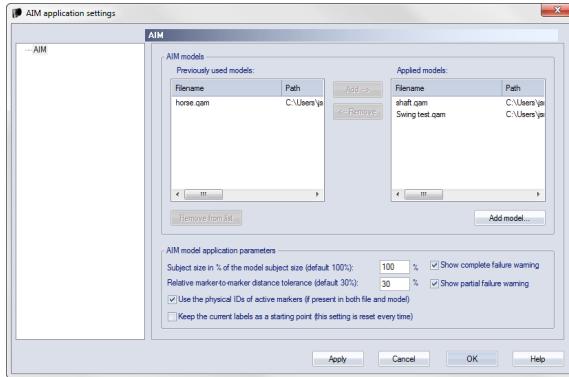

4. Check that the **Applied models** are correct. It is possible to apply several AIM models to the same file, for more information see “AIM models for multiple subjects” in the QTM user manual.
5. Click **OK**, the AIM module will try to apply the model to the measurement. Click **Cancel** in the dialog below to abort AIM.

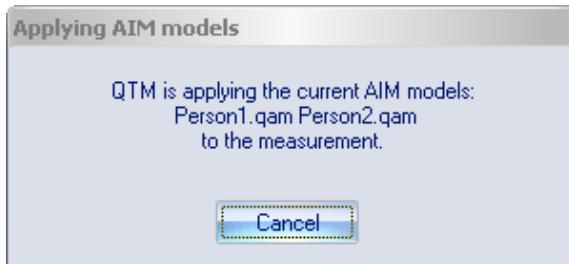

If any of the models cannot be applied to the trajectories the dialog below will appear. Showing how many of the bodies (models) were applied.

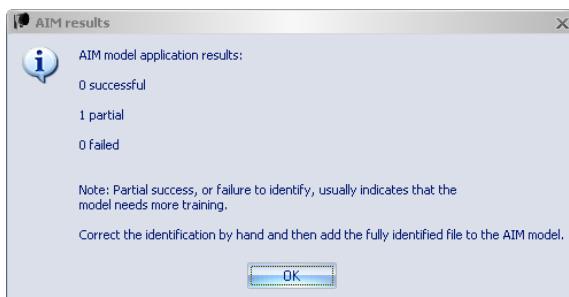

The **AIM results** dialog will display the result of all of the applied AIM models. The **Partial** results are AIM models where not all of the markers have been identified. For the **Failed** models none of the markers have been identified.

The most likely reason for the **Partial** result is that the model on the **AIM** page is incorrect for the captured motion. However, if the AIM model is correct it is recommended to manually identify the file and then add it to the existing AIM model, see “Generating an AIM model” on page 50.

It can also help to reduce the selected measurement range so that the AIM model is applied only on a smaller part of the measurement. For example, if the subject walks in and out of the volume it can help to reduce the selected measurement range to where the subject is inside the volume.

6. Check the file to see whether all of the trajectories have been identified correctly. There are two main reasons that trajectories are not identified, for more information about AIM see chapter "How AIM identifies trajectories" below.
 - AIM will join parts of trajectories that are within the AIM model definition and also split the trajectory if it is too much outside the AIM model definition. Therefore if there are some trajectories that are not identified in a part of the measurement it probably did not match the model in that part. The solution is to identify the markers and then add the file to the AIM model.
 - AIM sometimes cannot identify a trajectory if a nearby trajectory is missing during that part of the measurement.
7. When the AIM model has been successfully applied, the capture file must be saved to keep the changes.

When you apply a model as a processing step, either directly after a capture or in a batch process, it works exactly like when applying it manually. The model is set on the **AIM** page in the **Project options** dialog.

How AIM identifies trajectories

AIM identifies trajectories from angles and distances between markers. In the AIM model QTM has saved the angle and distance ranges of added measurement files. This is the reason why it works better when you add more measurements since the AIM model will then include more movements. When the model is applied AIM tries to find the solution that fits these ranges the best.

There are a few limitations to AIM that are good to know:

1. AIM does not work on only two markers. The main reason is that AIM works in local coordinates, i.e. it does not take into account whether "marker x is to the left of marker y". AIM only cares about "what location marker x has compared to the locations of markers y and z". This means that there is no left/right, over/below or in front of/behind in the model.
2. AIM needs to have movement in the first file used to create the model. Otherwise AIM will not identify the correct bones in the AIM model, since the bones are based on the shortest possible distances between markers.
3. AIM looks for the best solution in each frame and then joins the solutions. This means that AIM will join and split trajectories depending on the AIM model so that you get the best match. The measurement subject can leave and enter the volume and AIM will continue to join the trajectories.
4. Because AIM is working with relations between markers, all of the markers in the model file must be included in the measurement to ensure that all of the trajectories can be identified. E.g. if a marker on the shoulder disappears for a long time it is difficult for AIM to find the arm.

Qualisys AB is continuously working on improving AIM so please contact us if you have any questions.

Data export to other applications

The data in QTM can be exported to other applications for external analysis. The data export can be done either as a processing step or manually. To export the data manually go to the **File** menu and the **Export** sub-menu. The file formats that are

exported during processing are selected on the **Processing** page in the **Project options** dialog.

The data can be exported to the following formats:

TSV format

Tab separated file that can be opened in e.g. Microsoft Excel.

C3D format

Motion capture industry file format that can be used in e.g. Visual3D, Motion Monitor and MotionBuilder.

MAT file/directly to Matlab

The data is exported in a struct array either in a MAT-file or directly into Matlab.

DIFF format

A data file format primarily used in Japan.

IR radiation notice and environmental regulations

IR radiation notice

The Oqus camera uses short but quite strong infrared flashes to illuminate the markers. The flash is generated by LEDs on the front of the camera. The Oqus cameras are within the limits of the EN-60825-1-2001 Class 1 and the FDA CFR 1040.10 Class I classification, which means that the LED radiation is not considered to be hazardous.

However, any light of high intensity might be harmful to your eyes. Because infrared light is invisible to the human eye, you can be exposed to IR light without noticing. Therefore we recommend that you do not stare directly at the LEDs at a short distance for a prolonged time period when the camera is running.

EU customer information about about Waste Electrical and Electronic Equipment (WEEE)

In the European Union (EU), waste from electrical and electronic equipment (WEEE) is now subject to regulation designed to prevent the disposal of such waste and to encourage prior treatment measures to minimize the amount of waste ultimately disposed. In particular, the EU WEEE Directive 2002/96/EC requires that producers of electronic equipment be responsible for the collection, reuse, recycling and treatment of WEEE which the producer places on the EU market after August 13, 2005. Qualisys is providing the following collection process to comply with the WEEE Directive.

Qualisys WEEE Collection Process

If you have purchased Qualisys products in the EU on and after August 13, 2005, and are intending to discard these products at the end of their useful life, please do not dispose of them in a landfill or with household or municipal waste. Qualisys has labeled its electronic products with the WEEE label to alert our customers that products bearing this label should not be disposed with waste in the EU. Instead, Qualisys requests you to return those products using the instructions provided here, so that the products can be collected, dismantled for reuse and recycled, and properly disposed.

Qualisys will take back WEEE, i.e. all of the electrical equipment which is part of Qualisys equipment, from its customers within the EU. Please visit the website www.qualisys.com/weee or contact Qualisys AB at weee@qualisys.com for information on how to return your WEEE.

China ROHS (in Chinese)

有害物质声明

按照中华人民共和国电子工业标准SJ/T11364 2006的要求，此文档提供了由Qualisys AB生产制造的Oqus1-, 3-和5-系列的危险材料声明。

部件名称	有毒有害物质或元素					
	铅	汞	镉	六价铬	多溴联苯	多溴二苯醚
印制电路配件	x	o	o	o	o	o
显示器	x	o	o	o	o	o
按钮	o	o	o	o	o	o
内部配线	o	o	o	o	o	o
围栏	o	o	o	o	o	o
镜头	x	o	o	o	o	o
外接电缆及端口	o	o	o	o	o	o
交流电/直流电电源	o	o	o	o	o	o
纸质说明书	o	o	o	o	o	o

CD 说 明 书	o	o	o	o	o	o
-------------------	---	---	---	---	---	---

O: 表示该有毒有害物质在该部件所有均质材料中的含量均在SJ/T 11363-2006标准规定的限量要求以下。

X: 表示该有毒有害物质至少在该部件的某一均质材料中的含量超出SJ/T 11363-2006标准规定的限量要求。